

# Microarray Data Analysis Using BRB-ArrayTools Version 4.5.0 Beta1

Lori Long and Eric Polley  
May13-14, 2014

# Agenda

## Day 1

- I. What is BRB-ArrayTools?
- II. Installing BRB-ArrayTools and its required components
- III. Creating a collated project workbook
- IV. Data filtering and normalization options
- V. Scatterplot
- VI. Class Comparison
- VII. Gene Set comparison

# Agenda

## Day 2

- I. Clustering
- II. MDS
- III. Class Prediction
- IV. RNA-Seq Count Data Importer and DE analysis
- V. Plug-ins
- VI. Tutorial.
- VII. Hands-on.

Part I:

What is BRB-ArrayTools?

# BRB-ArrayTools

## An Integrated Software Tool for Microarray Data Analysis

- Developed under the direction of Dr. Richard Simon of the Biometric Research Branch, NCI.
- Software was developed with the purpose of deploying powerful statistical tools for use by biologists.
- Analyses are launched from user-friendly Excel interface. Also requires installation of R for running back-end programs. Current requirement for R is v 3.2.0. Publicly available from BRB website:

<http://linus.nci.nih.gov/BRB-ArrayTools.html>

# Features of BRB-ArrayTools

- Capability to collate (sort into an expression data matrix) microarray data from a set of experiments, and apply filtering and normalization.
- The focus of the software has been the implementation of statistical methodology which utilizes the sample descriptors (supervised analysis).
- Scatterplots, hierarchical clustering, and multidimensional scaling analyses also provide powerful visualization tools.
- Gene annotations are integrated into analysis output to inform the analysis results. Also, includes analyses using BioCarta, KEGG and Broad/MIT pathways.
- Advanced users may program their own plug-in analysis tools within BRB-ArrayTools.

# Limitations of BRB-ArrayTools

- Available only on the PC. As well as on an Apple macbook pro machine with Windows OS installed with Apple's bootcamp software/Parallels.
- Currently compatible with MS Vista/ Windows 7/Windows8 and Excel 2007/2010/2013.
- Importing and running analysis tools on large data sets may require a large memory capacity and be time-consuming.

# Recently Added new tools

- Dynamic Heatmap Viewer to generate interactive heatmaps on clustered data.
- Importing and Differential expression analysis on RNA-Seq count data.
- Importing Affymetrix ST gene chip data using the “oligo” package.



## Part II:

# Installing BRB-ArrayTools and its required components

# Installing BRB-ArrayTools

- <http://linus.nci.nih.gov/BRB-ArrayTools.html>
- Register to obtain a user name and password by going to the guestbook.
- Select the version you wish to download.
- Currently available BRB-ArrayToolsv4.4.0.
- BRB-ArrayTools v4.5.0 Beta1 is coming soon.

# BRB-ArrayTools Web Page

## BRB-ArrayTools

**Developed by: Richard Simon & BRB-ArrayTools Development Team**

BRB-ArrayTools is an integrated package for the visualization and statistical analysis of DNA microarray gene expression data. It was developed by professional statisticians experienced in the analysis of microarray data and involved in the development of improved methods for the design and analysis of microarray based experiments. The array tools package utilizes an Excel front end. Scientists are familiar with Excel and utilizing Excel as the front end makes the system portable and not tied to any database. The input data is assumed to be in the form of Excel spreadsheets describing the expression values and a spreadsheet providing user-specified phenotypes for the samples arrayed. The analytic and visualization tools are integrated into Excel as an add-in. The analytic and visualization tools themselves are developed in the powerful R statistical system, in C and Fortran programs and in Java applications. Visual Basic for Applications is the glue that integrates the components and hides the complexity of the analytic methods from the user. The system incorporates a variety of powerful analytic and visualization tools developed specifically for microarray data analysis.



[BRB-ArrayTools:  
Gene expression data analysis](#)




[FAQs & Answers](#)  
[BRB-ArrayTools Message Board](#)



[BRB-ArrayTools Data Archive  
for Human Cancer Gene Expression](#)



[Email BRB-ArrayTools Support](#)

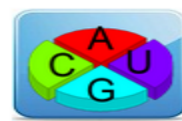

[BRB Digital Gene Expression:  
Pre-Process RNA-Seq Files](#)



[Book for DNA Microarray Analysis](#)



[Publications Based on BRB-ArrayTools  
Analyses](#)



[BRB-ArrayTools User Community  
Institution List](#)

[Guest Book](#) | [Message Board](#) | [Licenses](#) | [Reprints and Presentations](#)

# Full Installation File

- This file is a bundle of all the necessary components. Rv3.2.0 and java are included along with ArrayTools and CGHTools.



## Biometric Research Branch

*Division of Cancer Treatment and Diagnosis*

### BRB-ArrayTools Licensing Agreement

Licensing agreements for BRB-ArrayTools differ for U.S. Government users, academic/non-profit users, or commercial users. All users must agree to the following conditions:

1. All publications based on BRB-ArrayTools analyses will contain the acknowledgment: "Analyses were performed using BRB-ArrayTools developed by Dr. Richard Simon and BRB-ArrayTools Development Team."
2. The package or any of its components will not be distributed to others without prior approval of Dr. Richard Simon.
3. The package is offered without support and the user will not hold the National Cancer Institute, SRA International or the EMMES Corporation liable for any damages resulting from the use of BRB-ArrayTools.

**If you are a new user, please complete registration at our [Guest Book](#) before installation. This software is free for non-commercial use. Commercial users should contact Michael Shmilovich at [shmilovichm@od.nih.gov](mailto:shmilovichm@od.nih.gov) or (301)435-5019 for licensing. Commercial users can use a 60-day trial version by registering at our guest book.**



### [Download BRB-ArrayTools version 4.4.0 Stable release](#)

(released on Nov. 20, 2014)

- [What's New](#)
- [Instructions for Excel 2007/2010 Users](#) to set security level and load the Add-Ins into Excel 2007/2010 after installation
- [New Scatter Plot demo video](#)


The following documentation files are included in the above software installations, or may be downloaded separately for perusal prior to installation of the software.

- [Readme file](#)
- [Manual file](#)
- [Galaxy instruction file](#)

# Installing BRB-ArrayTools


- The ALL in ONE file called “ArrayTools\_v4\_5\_0\_Beta\_1.exe”.
- Includes installation of ArrayTools, R 3.2.0, Java, RServe package, Microsoft Visual C++ sp1 redistribution package.

# Installing BRB-ArrayTools

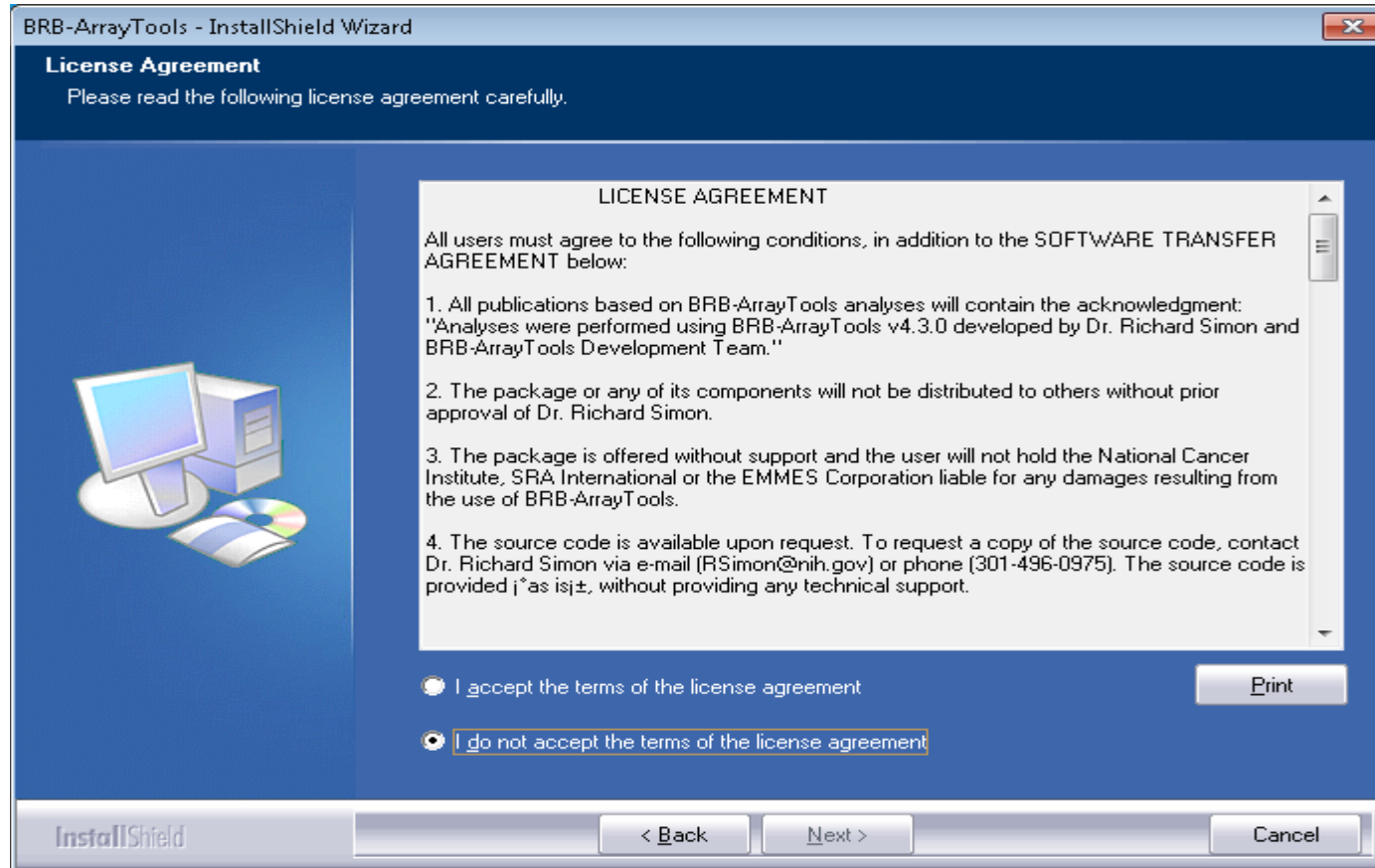


# Installing BRB-ArrayTools

- Select installation folder.

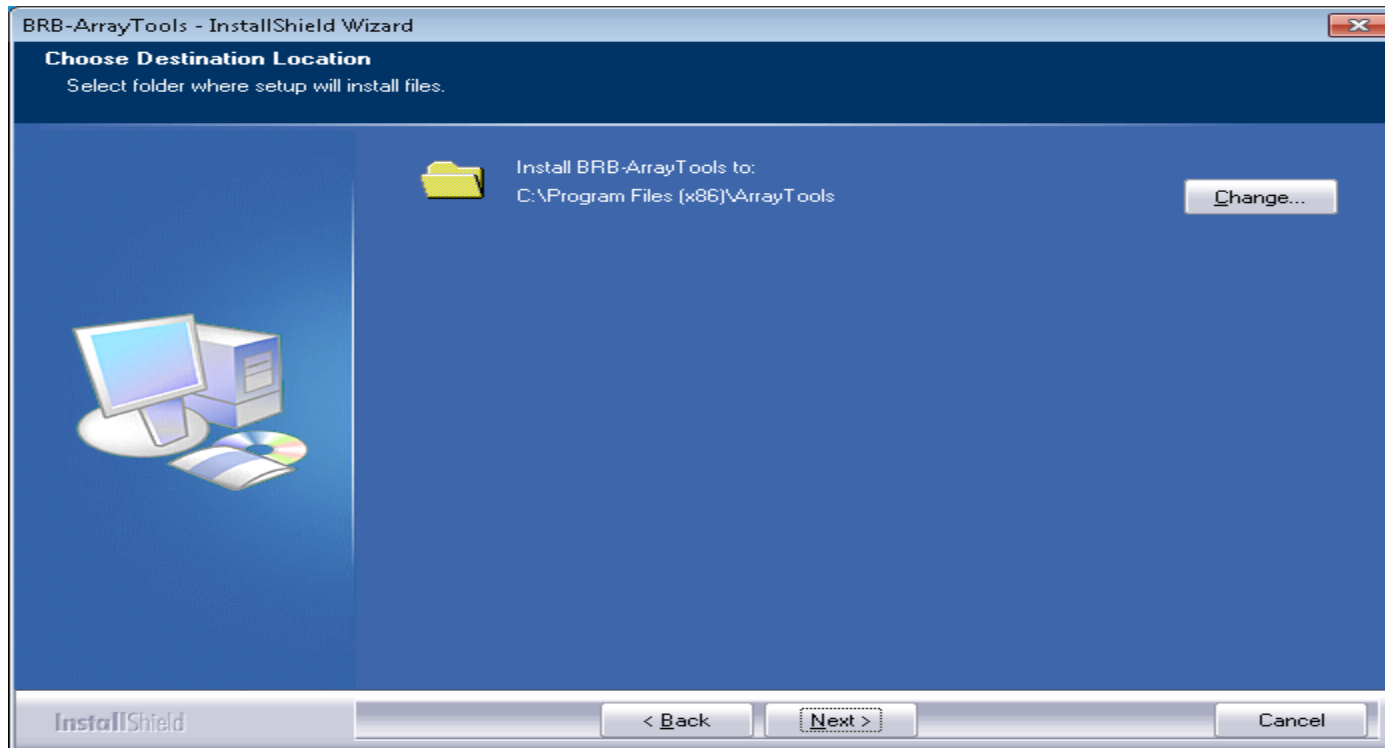

# Installing BRB-ArrayTools

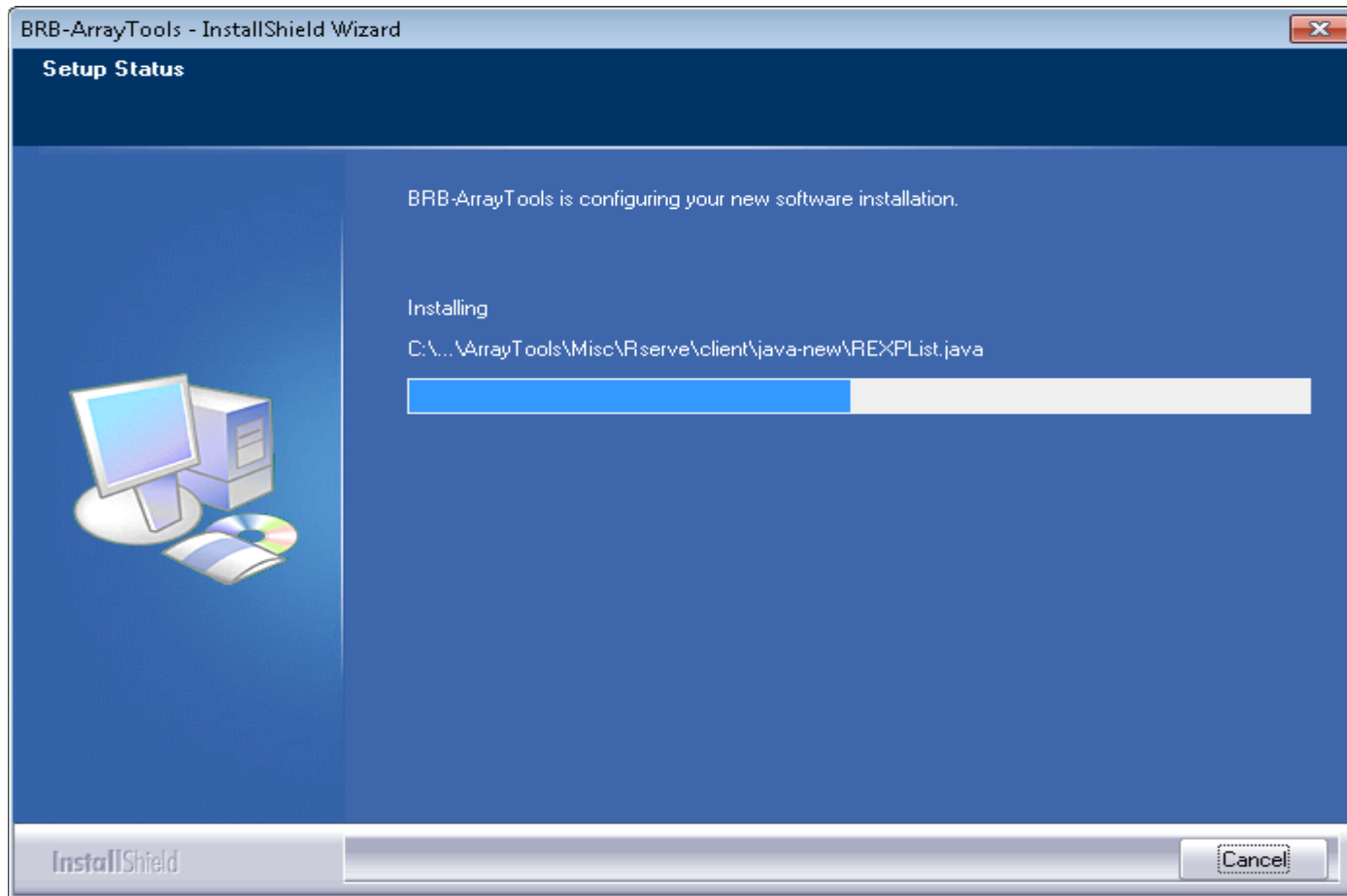




# Installing BRB-ArrayTools

- Proceed to install required components



# Installing BRB-ArrayTools

- After successfully installing BRB-ArrayTools, you will be prompted with the message below.
- Click “OK” as the software has been installed as an add-in to Excel.

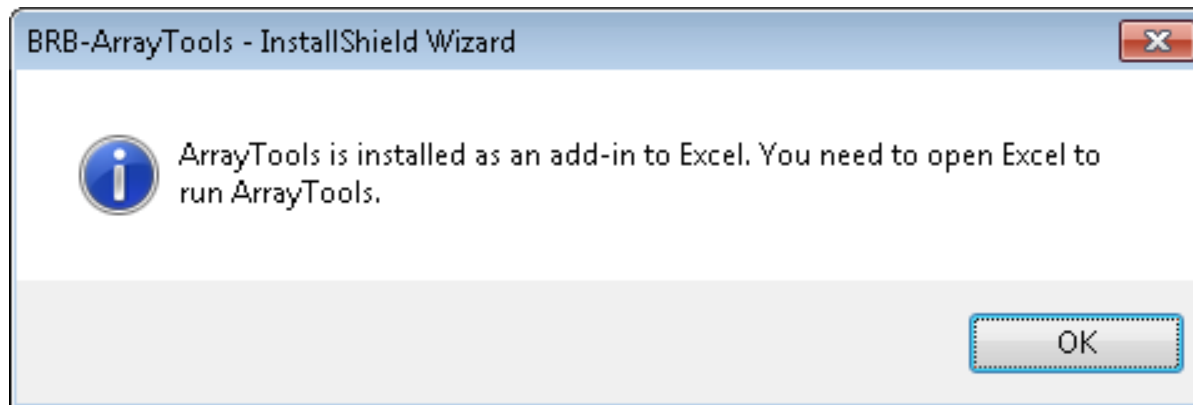

# Excel 2007/2010- loading the add-in

- 1: Click on the Microsoft 'Office' button on the top left corner of the Excel menu.
- 2. Then, select the "Excel Options" button on the bottom right.
- 3: Click on "Trust Center"
- 4. Then click on "Trust Center Settings"
- 5: Choose the "Macro Settings" from the left hand panel.
- 6. Check "Enable all macros" and "Trust access to VBA project."
- 7. Click the "OK" button.
- 8: Choose the "Add-ins" option from the left hand tab.
- 9. Click "BRB-ArrayTools" on the Active or Inactive application add-in.
- 10. Hit the "Go" button down at the bottom.
- 11. Check all the three "Add-ins", BRB-ArrayTools and CGHTools.
- 12. Then click OK.
- If you don't see the "Add-ins" ribbon along side "Home Insert...Review View" panel at the top then please close Excel and re-start.
- On clicking on Add-Ins tab, the Add-Ins should be listed there namely: ArrayTools and CGHTools add-ins.

# [Hands-on instructions]

## [Getting started]


1. Open **Excel**.
2. Click on **Tools → Add-ins**, and see that **BRB-ArrayTools** is loaded as an add-in.
3. When BRB-ArrayTools is loaded as an add-in, you will find an **ArrayTools** menu. This is the interface for all BRB-ArrayTools functions.
4. Click on **ArrayTools → Getting started**.
5. Here you will see the **Tutorial** and **Open a sample dataset** options.
6. For Office 2007, click on the “Add-ins” and you should find “ArrayTools”.


# ArrayTools Menu

- From the “ArrayTools” pull down menu, you can navigate to the following:
  - Manuals
  - Tutorials
  - Utilities
  - Support
  - License and version information

## Part III:

Getting your data into  
BRB-ArrayTools:  
Creating a project workbook

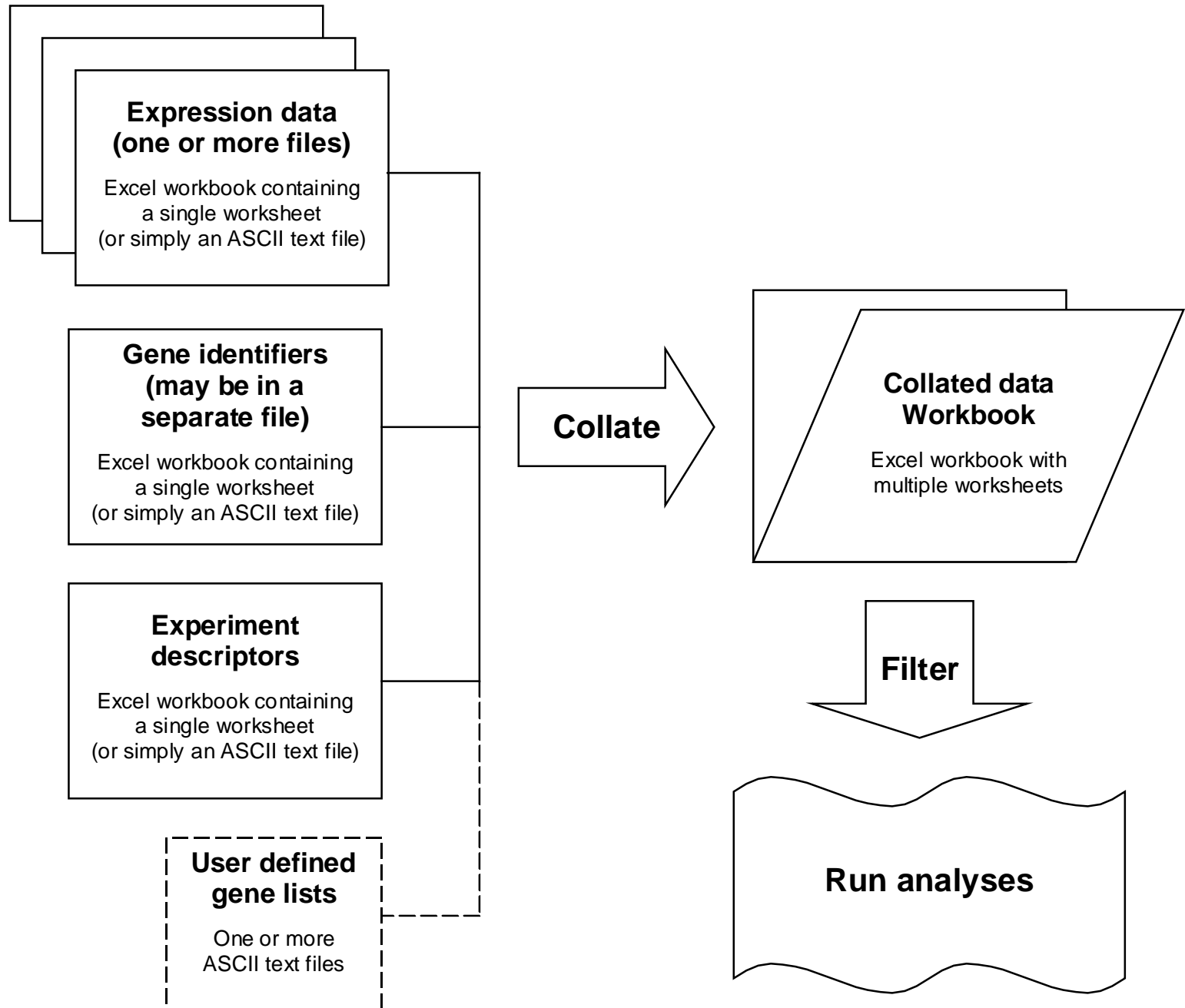


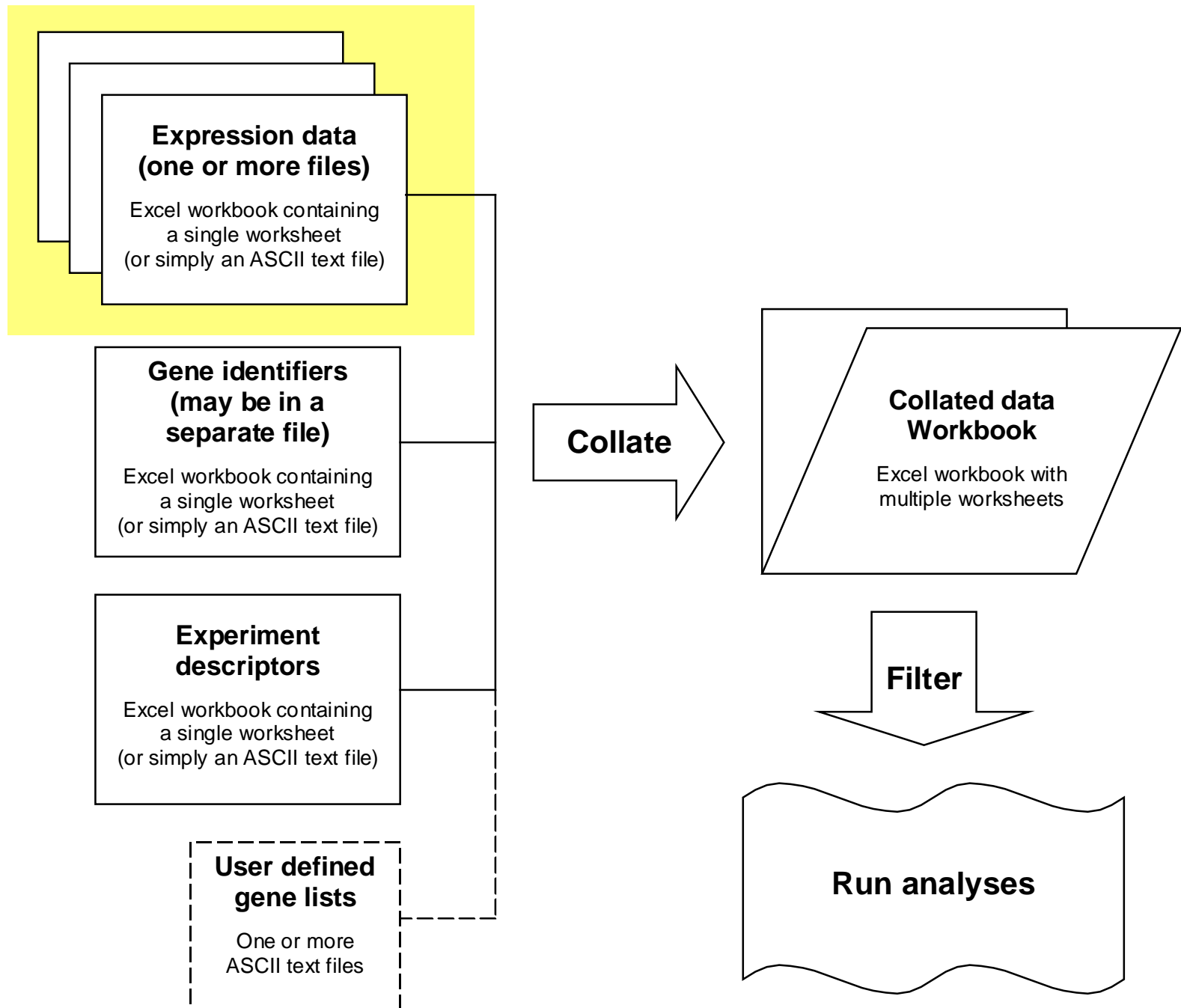



# Expression data

- Input data as tab-delimited ASCII files (or Excel spreadsheets) in one of the following formats:
  1. Horizontally aligned
  2. Separate files
- Files may contain expression data in the form of signal (or single-channel expression summary), dual-channel intensities, or expression ratios (for dual-channel data). Data may or may not have been already log-transformed. Flags, detection call, and spot size may also be used. All other variables will be ignored.

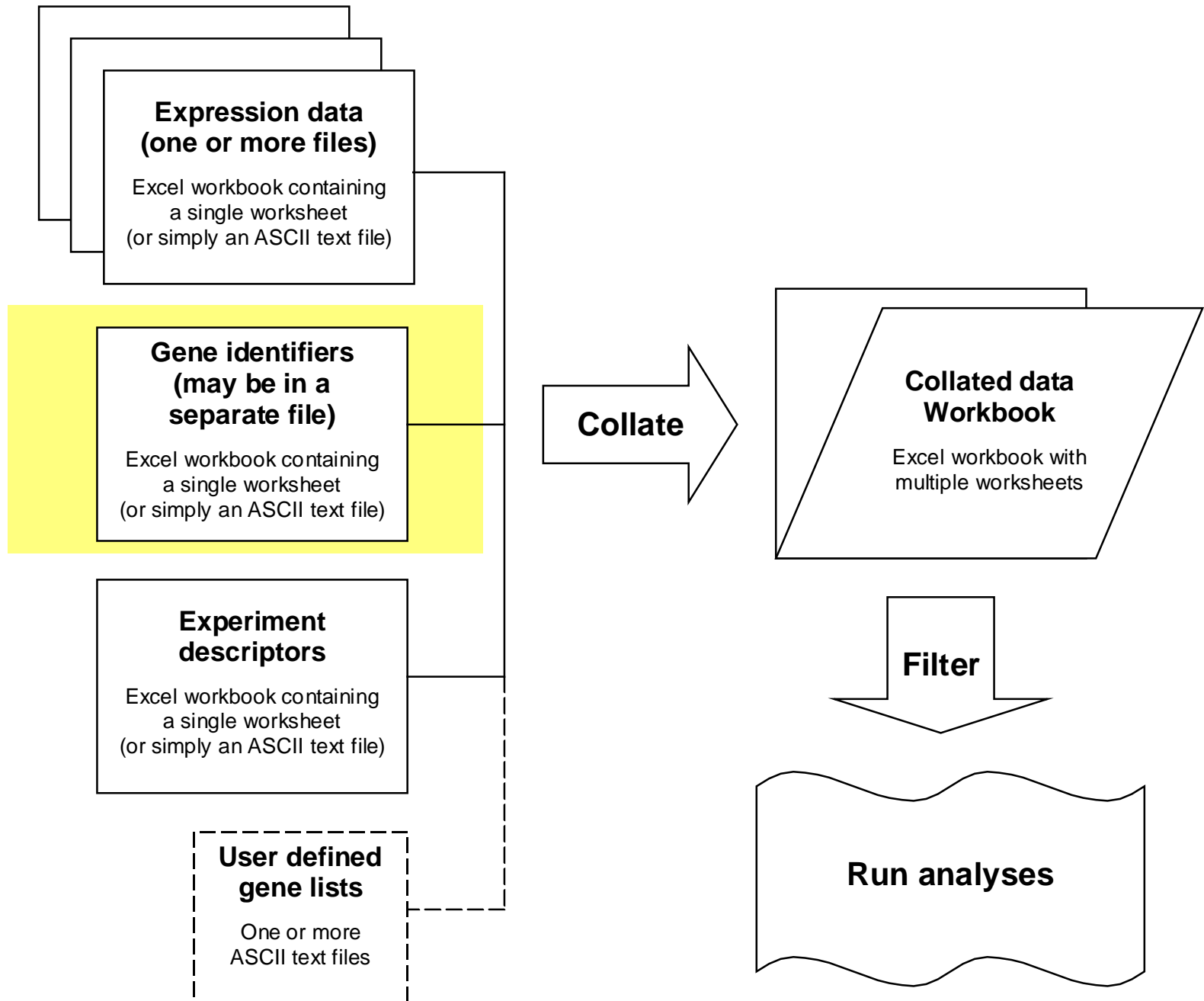

# Gene identifiers


- A gene identifiers file is optional, but highly recommended for annotation purposes.
- Gene identifiers which may be imported are: clone IDs, UniGene cluster IDs, gene symbols, GenBank/RefSeq accessions, Entrez IDs, Ensembl IDs and probe set IDs (for Affymetrix arrays). For microRNA data, the microRNA IDs can also be imported.

# Gene identifiers

Two examples of a gene identifier file

	A	B	C	D	E	F	G
1	UniqueID	Ensembl	Name	Symbol	EntrezID	Accession	UGCluster
7	ENSG00000000938	ENSG00000000938	FGR proto	FGR	2268	NM_00104	Hs.1422
8	ENSG00000000971	ENSG00000000971	compleme	CFH	3075	NM_00018	Hs.363396
19	ENSG00000001630	ENSG00000001630	cytochrom	CYP51A1	1595	NM_00078	Hs.417077
23	ENSG00000002330	ENSG00000002330	BCL2-asso	BAD	572	NM_00432	Hs.370254
26	ENSG00000002587	ENSG00000002587	heparan s	HS3ST1	9957	NM_00511	Hs.507348
28	ENSG00000002745	ENSG00000002745	wingless-1	WNT16	51384	NM_01608	Hs.272375
29	ENSG00000002746	ENSG00000002746	HECT, C2 a	HECW1	23072	NM_00128	Hs.164453
36	ENSG00000003137	ENSG00000003137	cytochrom	CYP26B1	56603	NM_00127	Hs.91546
42	ENSG00000003436	ENSG00000003436	tissue fact	TFPI	7035	NM_00103	Hs.516578
52	ENSG00000004468	ENSG00000004468	CD38 mole	CD38	952	NM_00177	Hs.479214
53	ENSG00000004478	ENSG00000004478	FK506 bin	FKBP4	2288	NM_00201	Hs.524183
56	ENSG00000004660	ENSG00000004660	calcium/c	CAMKK1	84254	NM_03225	Hs.8417

	ProbeSet	Name	Accession	UGCluster	Symbol	EntrezID	Chromosome	Cytoband
1	16966685	spermato	NM_14526	Hs.527090	SPATA18	132671	4	4q12
2	16966855	v-kit Hard	NM_00022	Hs.479754	KIT	3815	4	4q12
3	16967376	casein alp	NM_00185	Hs.3155	CSN1S1	1446	4	4q21.1
4	16967631	solute car	NM_00105	Hs.5462	SLC4A4	8671	4	4q21
5	16968356	bone mor	NM_00120	Hs.387411	BMP3	651	4	4q21
6	16968447	1-acylglyc	NM_03271	Hs.99196	AGPAT9	84803	4	4q21.23
7	16968675	heat shock	AY956764	Hs.744112	HSP90AB3	3327	4	4q22.1
8	16968680	secreted p	NM_00125	Hs.313	SPP1	6696	4	4q22.1
9	16968878	multimeri	NM_00735	Hs.268107	MMRN1	22915	4	4q22
10	16969229	dual adap	NM_01435	Hs.436271	DAPP1	27071	4	4q25-q27

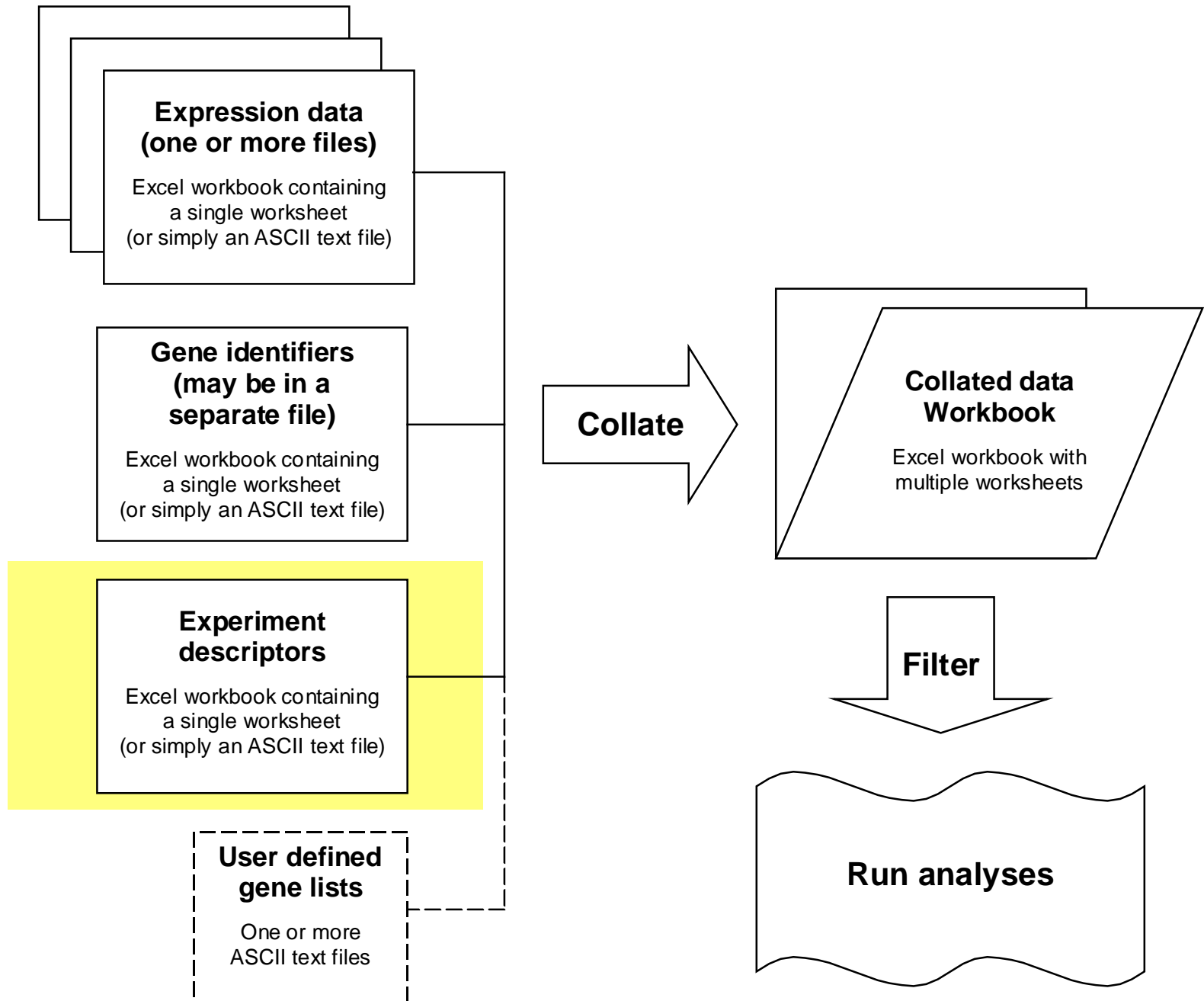


# Experiment (Array) descriptors

- An experiment descriptors file describes the samples used for each array, and is mandatory.
- After the header row, each row in this file represents one array or sample, and each column represents one descriptor variable.
- First column contains array id, which is matched against file names when expression data is in separate files format.
- Subsequent columns contain descriptions, phenotype class labels, patient outcome, and other sample or experiment information.
- The descriptor variable columns may include information such as: patient ids, class labels, technical replicate indicators, reverse fluor indicators, and other variables used for labeling purposes.
- A COPY of the original experiment descriptor file will appear in the experiment descriptor sheet of the collated project workbook. The experiment descriptor sheet in the collated project workbook may be further edited as you analyze the data.

# Experiment descriptors

## Describes the samples used for each array

A	B	C	D	E	F	G	H	I	J	K	L	M	N	O	P	Q	R	S	T	U
Experiment Name	Type	Annotation	Full annot:	Clinical P	Path II	Max tumor	Tumor per	Tumor per	Operation	Age	Survival	Status at 1	Status	Sex	Smoking	Stage AJC	Summary	Site of rel	Censor	
AD-043T2-A7-1	AD	AD043	043T2-A7-	ad m (met adm		40	40		bilobectomy	47.58904	14.1	1 d		M	2.5 M1					Uncensored
AD-111T2-A8-1	AD	AD111	111T2-A8-	ad 2.0 m-p		80	80		wedge	76.03562	72.4	0	1 F		40 T1Nxm	IA				Censored
AD-114T1-A9-1	AD	AD114	114T1-A9-	ad 1.5 m-w w/BAC		35	5	35	wedge	73.47945	66.2	0	1 M		40 T1NOMx	IA				Censored
AD-115T1-A12-1	AD	AD115	115T1-A12-	ad 6.5 m adm/adw		40	40	15	lobectomy	70.45479	21.9	1	3 F		75 T2N1M0	IIB	lung, LN			Uncensored
AD-118T1-A13-1	AD	AD118	118T1-A13-	ad 2.5 m adm		60	60	60	lobectomy	68.58356	49.6	1	3 M		25 T1NOMx	IA	lung, LN			Uncensored
AD-119T3-A195-8	AD	AD119	119T3-A19-	ad 3.5 w Eadp		50	25	50	lobectomy	58.79452	76.6	0	1 F		10 T2NOMx	IB				Censored
AD-120T1-A226-8	AD	AD120	120T1-A22-	ad 8.0 m adm		65	30	65	lobectomy	68.35616	38.9	1	3 M		54 T2NOMx	IB	bone			Uncensored
AD-122T3-A197-8	AD	AD122	122T3-A19-	ad 5.0 m adm		65	50	65	lobectomy	72.65753	33.9	1	3 F		0 T2N1Mx	IIB	lung			Uncensored
AD-123T1-A25-1	AD	AD123	123T1-A25-	ad 5.0 m adm,pap		60	50	60	lobectomy	60.35616	74	0	1 M		126 T3NOMx	IIB				Censored
AD-127T1-A14-1	AD	AD127	127T1-A14-	ad 1.8 p adp		80	70	80	lobectomy	65.49041	8.2	1	3 F		69 T1N2Mx	IIIA	LN			Uncensored
AD-130T1-A1-1	AD	AD130	130T1-A1-	ad 15.0 B BAC		80	50	80	lobectomy	74.71233	7.1	1 d		M	100 T2N1Mx	IIB				Uncensored
AD-131T1-A15-A200-1	AD	AD131	131T1-A15-	nsc 2.1 p adp		80	80	80	wedge	59.63562	76	0	1 F		100 T1NOMx	IA				Censored
AD-136T2-A201-8	AD	AD136	136T2-A20-	ad 4.0 m adm		75	70	75	lobectomy	65.6	31.4	0	1 F		100 T2NOMx	IB				Censored
AD-157T1-A246-10	AD	AD157	157T1-A24-	ad poor (r adp	poor, ac	70	50	70	pleurectomy	73.15342	4.1	1 d		F	55 M1					Uncensored
AD-158T1-A247-10	AD	AD158	158T1-A24-	ad 3.9 m adw,acina		33	33	20	lobectomy	64.00548	40.7	1	3 M		107.5 T2NOMx	IIB	lung, bone			Uncensored
AD-159T1-A229-8	AD	AD159	159T1-A22-	ad 5.5 m adw,acina		60	55	60	lobectomy	70.50685	19.7	1 d		M	80 T2N1Mx	IIB				Uncensored
AD-162T2-A230-8	AD	AD162	162T2-A23-	ad 3.5 m admod,aci		60	60	60	lobectomy	74.87123	41.7	0	1 F		60 T2NOMx	IB				Censored
AD-163T1-A203-8	AD	AD163	163T1-A20-	ad 3.2 m admod,aci		80	70	80	lobectomy	69.55616	86.7	1 d		F	100 T2N1Mx	IIB				Uncensored
AD-164T1a-A206-8	AD	AD164	164T1a-A2-	ad 4.5 p adp	poor, ac	90	90	70	lobectomy	68.15342	15	1	3 M		80 T3NOMx	IIB	LN			Uncensored
AD-167T1 (d dup)-A210-8	AD	AD167	167T1-A21-	ad 2.5 w w adw,acina		75	70	75	lobectomy	76.64384	41.7	0	1 M		0 T2NOMx	IB				Censored
AD-169T2-A211-8	AD	AD169	169T2-A21-	ad 2.5 m adw/pap o		80	70	80	lobectomy	46.97808	20	1	3 F		21.6 T2NOMx	IB	bone, myc			Uncensored
AD-170T1-A251-10	AD	AD170	170T1-A25-	ad 2.5 w w BAC & pa		85	35	85	lobectomy	60.94521	78.4	0	1 F		60 T1NOM0	IA				Censored
AD-172T2-A213-8	AD	AD172	172T2-A21-	ad 3.0 m Eadmod,aci		70	70	70	lobectomy	68.90137	66.7	0	1 F		35 T2N1Mx	IIB				Censored
AD-173T1a-A23-1	AD	AD173	173T1a-A2-	ad 5.0 m admod,aci		70	65	70	lobectomy	56.78356	22.3	1 d		F	27 T2N1Mx	IIB				Uncensored
AD-177T1-A21-1	AD	AD177	177T1-A21-	ad 3.2 p admod,aci		50	0	50	bilobectomy	59.44658	2.8	1	4 M		100 T3N1Mx	IIIA				Uncensored
AD-178T2-A22-1	AD	AD178	178T2-A22-	ad 1.9 w Eadmod,aci		80	80	80	lobectomy	68.94795	29	1	3 M		80 T1NOMx	IA	lung, LN, t			Uncensored

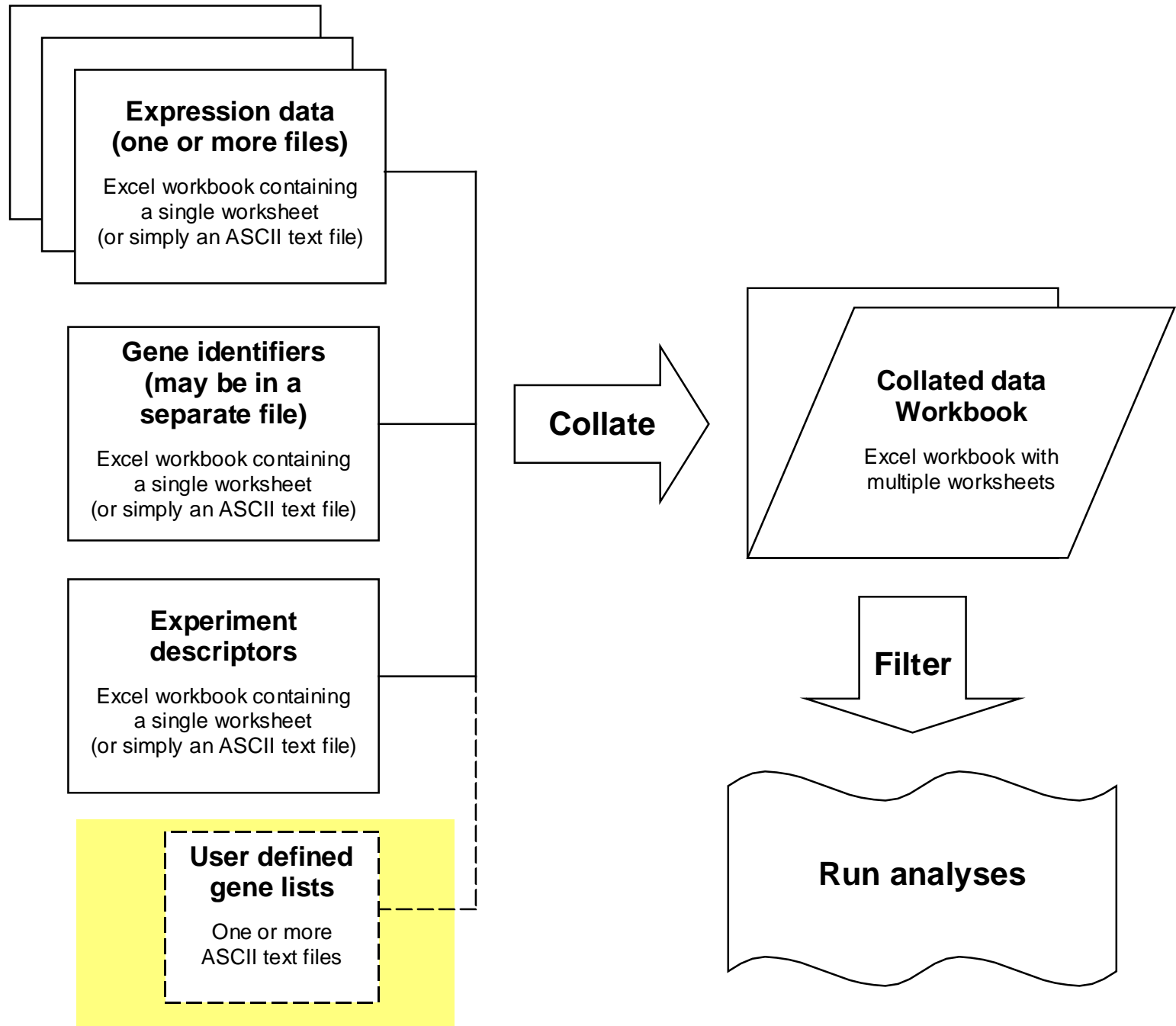





# Gene lists

- Genelists are used for annotation and for defining subsets for analysis. These files are located in the ArrayTools installation folder.
- Two types of genelists: Preloaded with BRB-ArrayTools, and user-defined
- User-defined genelists are simply text files which the user creates, containing a label specifying the type of identifier, followed by a list of gene identifiers. The file should be appropriately named to indicate what type of genes are in the list. Some user-defined genelists are automatically produced as the result of an analysis, such as class comparison, class prediction, survival analysis, and hierarchical clustering of genes.
- User-defined genelists are stored in the “project” folder (for project specific) or ArrayTools folder (visible to all projects.)

# Genelists pre-loaded with BRB-ArrayTools

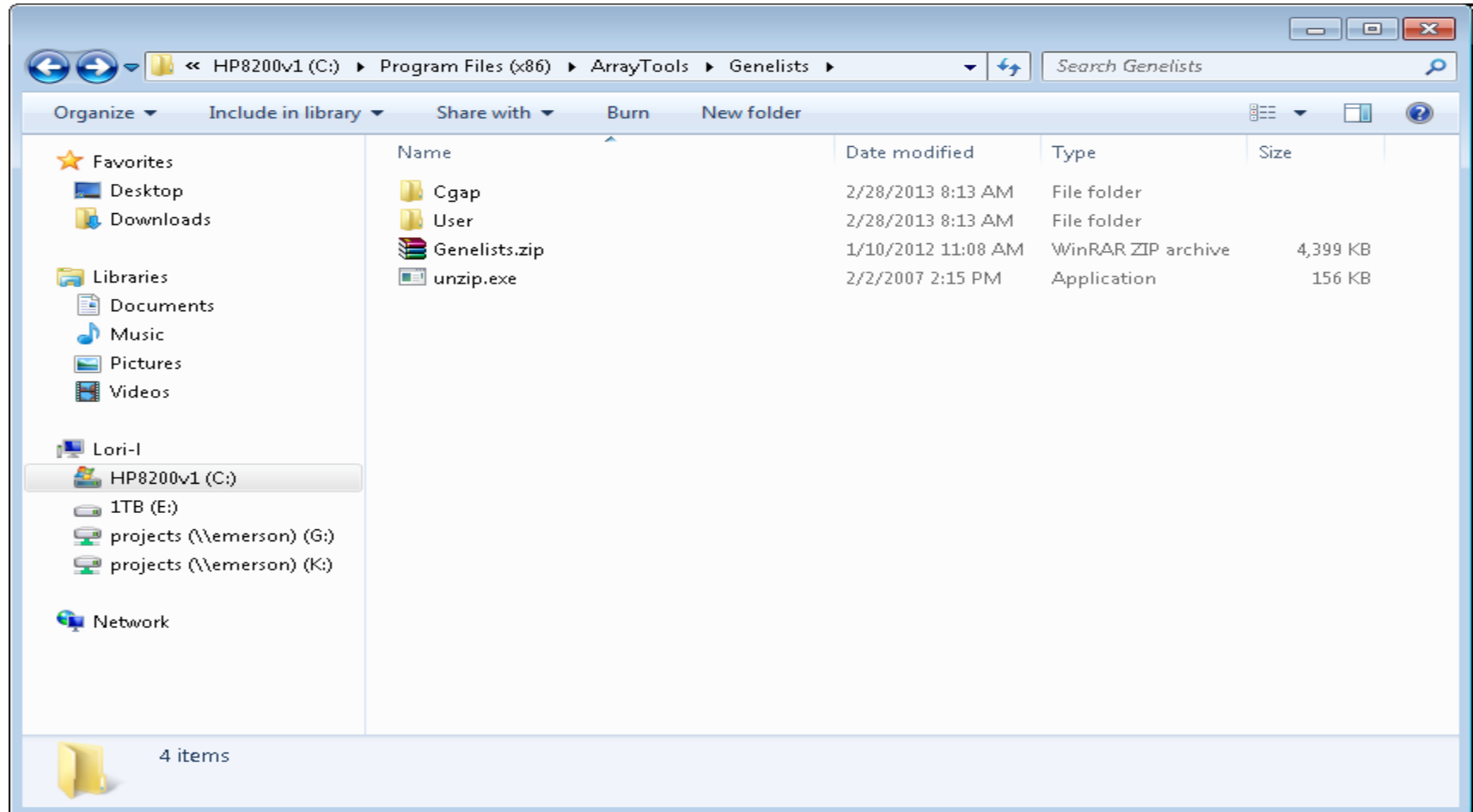

# Genelists Folder

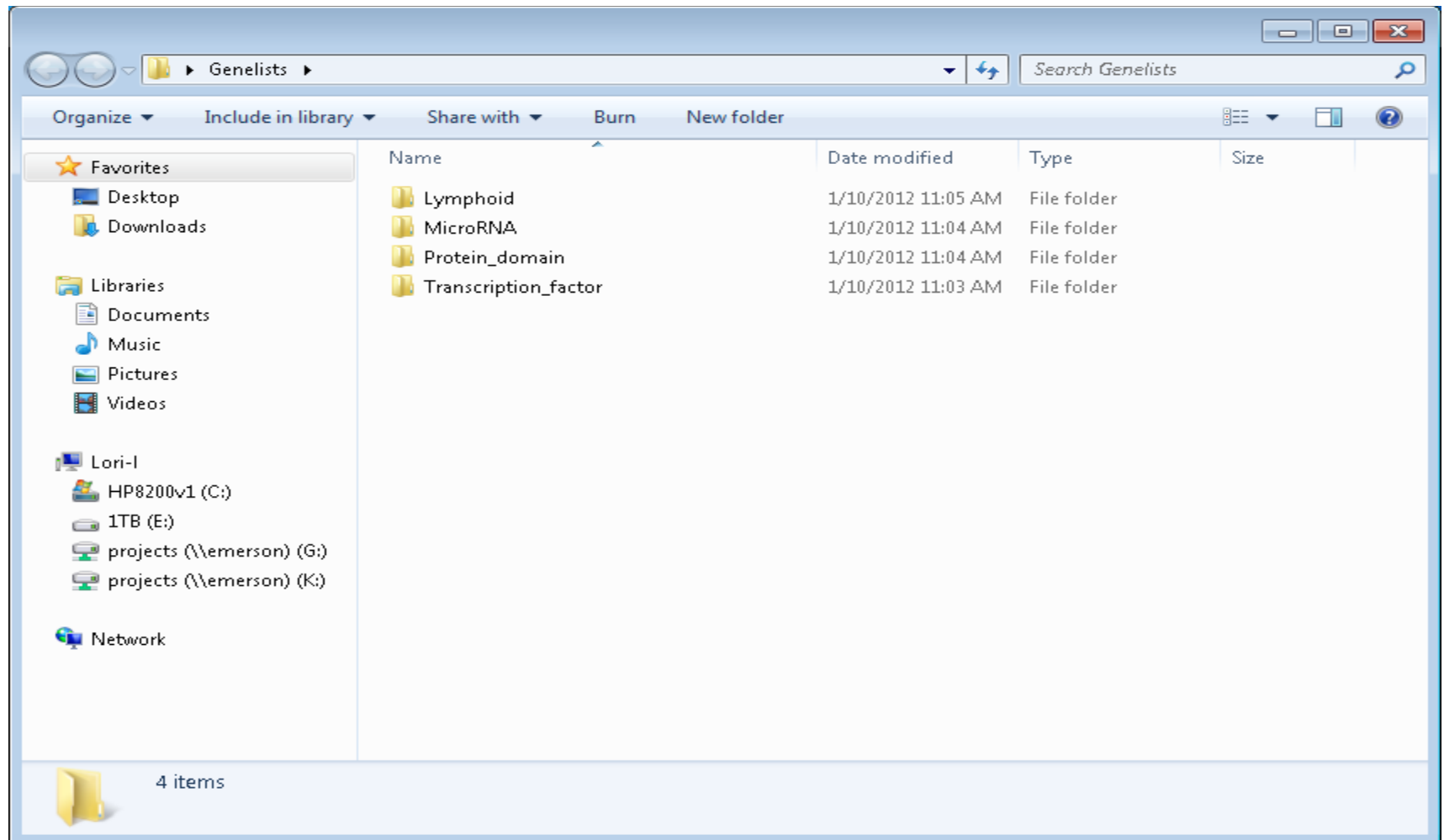

# Gene lists

## Cancer Genome Anatomy Project

CGI: Angiogenesis - Microsoft Internet Explorer

Address <C:\Program Files\ArrayTools\Genelists\CGAP\angiogenesis.html> Go

### Angiogenesis

- This collection curated by Elise Kohn ([ek1b@nih.gov](mailto:ek1b@nih.gov))

Gene	Description	Sequences	Sequence assembly	Predicted SNPs having score $\geq 0.99$
<a href="#">ADM</a>	Adrenomedullin	D14874	<a href="#">D14874</a>	<a href="#">1</a>
<a href="#">ANG</a>	Angiogenin, ribonuclease, RNase A family, 5	M11567	<a href="#">M11567</a>	<a href="#">1</a>
<a href="#">ANGPT1</a>	Angiopoietin 1	D13628, U83508	<a href="#">AF004327</a>	
<a href="#">ANGPT2</a>	Angiopoietin 2	AF004327		
<a href="#">ANGPT3</a>	Angiopoietin 3	AF107253		
<a href="#">ANGPT4</a>	Angiopoietin 4	AF113708		
<a href="#">ANPEP</a>	Alanyl (membrane) aminopeptidase (aminopeptidase N, aminopeptidase M, microsomal aminopeptidase, CD13, p150)	M22324	<a href="#">M22324</a>	<a href="#">2</a>
<a href="#">ARNT</a>	Aryl hydrocarbon receptor nuclear translocator	M69238		
<a href="#">BDK</a>	Bradykinin			
<a href="#">BDKRB2</a>	Bradykinin receptor B2	M88714, X86162, X86172, X86173	<a href="#">X86163</a>	<a href="#">1</a>


My Computer

# Gene lists

## User-defined text files

ClassComparison.txt - Notepad

ProbeSet	Accession	UGCluster	Symbol	EntrezID
S82024_at	S82024	Hs.521651	STMN2	11075
U52828_s_at	U52828	Hs.314543	CTNND2	1501
D21267_at	D21267	Hs.167317	SNAP25	6616
M14483_rnal_s_at	M14483	Hs.459927	PTMA	5757
U38810_at	U38810	Hs.584776	MAB21L1	4081
D82347_at	D82347	Hs.574626	NEUROD1	4760
M93119_at	M93119	Hs.89584	INSML	3642
X86809_at	X86809	Hs.517216	PEA15	8682
L10373_at	L10373	Hs.441664	TSPAN7	7102
U50822_rnal_s_at	U50822	Hs.574626	NEUROD1	4760
U76421_at	U76421	Hs.474018	ADARB1	104
U29195_at	U29195	Hs.3281	NPTX2	4885
U00802_s_at	U00802	Hs.130316	DBN1	1627
U96136_at	U96136	Hs.314543	CTNND2	1501
M33308_at	M33308	Hs.643896	VCL	7414
X02761_s_at	X02761	Hs.203717	FN1	2335
HG3454-HT3647_at	HG3454-HT3647			
D49958_at	D49958	Hs.75819	GPM6A	2823
S45630_at	S45630	Hs.53454	CRYAB	1410
M25667_at	M25667	Hs.134974	GAP43	2596
L48513_at	L48513	Hs.530077	PON2	5445
M29536_at	M29536	Hs.429180	EIF2S2	8894
M97287_at	M97287	Hs.517717	SATB1	6304
X70476_at	X70476	Hs.75724	COPB2	9276
U25034_s_at	U25034	Hs.504703	NNAT	4826
L10338_s_at	L10338	Hs.436646	SCN1B	6324
X51405_at	X51405	Hs.75360	CPE	1363
S78296_at	S78296	Hs.500916	INA	9118
M98539_at	M98539	Hs.446429	PTGDS	5730
U30521_at	U30521	Hs.36053	C5orf13	9315
M37457_at	M37457	Hs.515427	ATPLA3	478
AJ001421_at	AJ001421	Hs.525527	RER1	11079
X05196_at	X05196	Hs.155247	ALDOC	230
HG658-HT658_f_at	HG658-HT658			
HG3236-HT3413_f_at	HG3236-HT3413			
U45955_at	U45955	Hs.495710	GPM6B	2824
M94250_at	M94250	Hs.82045	MDK	4192
M93426_at	M93426	Hs.489824	PTPRZ1	5803
HG3597-HT3800_f_at	HG3597-HT3800			
U90915_at	U90915	Hs.433419	COX4I1	1327
U60975_at	U60975	Hs.368592	SORL1	6653
Y09836_at	Y09836	Hs.335079	MAP1B	4131
S54005_s_at	S54005	Hs.446574	TMSB10	9168
D12676_at	D12676	Hs.349656	SCARB2	950
D78012_at	D78012	Hs.135270	CRMP1	1400
U48705_rnal_s_at	U48705	Hs.631988	DDR1	780
U76638_at	U76638	Hs.591642	BARD1	580
J04173_at	J04173	Hs.632918	PGAM1	5223

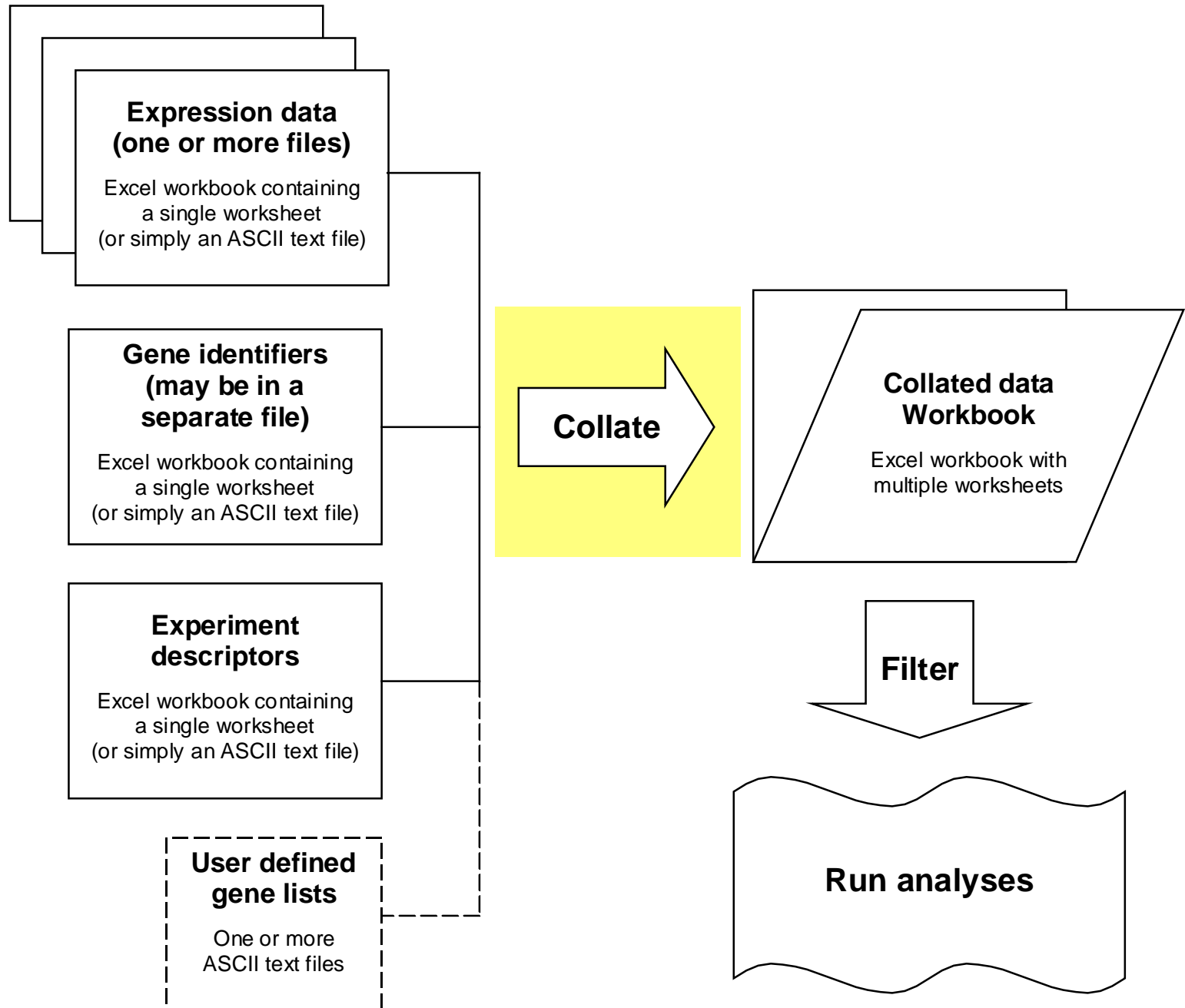

# Data importers

- General format data: The general format importer can be used to import most data formats.
- Data Import Wizard: This importer can be used to import a variety of data types.
- NCBI GEO GDS data Importer: This importer can be used to import GDS datasets from GEO archive.
- NCBI GEO GSE Data Importer.
- Affymetrix gene 1.0 st array importer.
- RNA-Seq count data importer.

# Data types

- The following data types can be imported through Data Import Wizard:
  - 3'-IVT array data in .CEL file format
  - Affymetrix gene 1.0 ST array data in .CEL file format
  - Affymetrix probe-set summary level data
  - Illumina expression data
  - Illumina methylation data
  - Agilent single and dual-channel data
  - Genepix single and dual-channel data
  - mAdb archive data
  - RNA-Seq data pre-processed by Galaxy




# Pomeroy Dataset

- **On the Desktop**, browse for the folder called “**Pomeroy**”.
- In this folder there are two files namely:  
**Dataset\_A2\_multiple\_tumor\_samples.txt**  
**ExpDescrMedulo.xls**
- The **Dataset\_A2\_multiple\_tumor samples.txt** contains the raw expression MAS5.0 summary values for all the arrays.
- The **ExpDescrMedulo.xls** contains the experiment descriptor file.

# [Hands-on instructions]

## [Importing Pomeroy Data set]

- Click on **ArrayTools** → **Getting started** → **Data Import Wizard**
- Select the option from the pull down menu- “**Affymetrix probeset-summary data**”.
- Choose the option that the expression data is combined into one file.

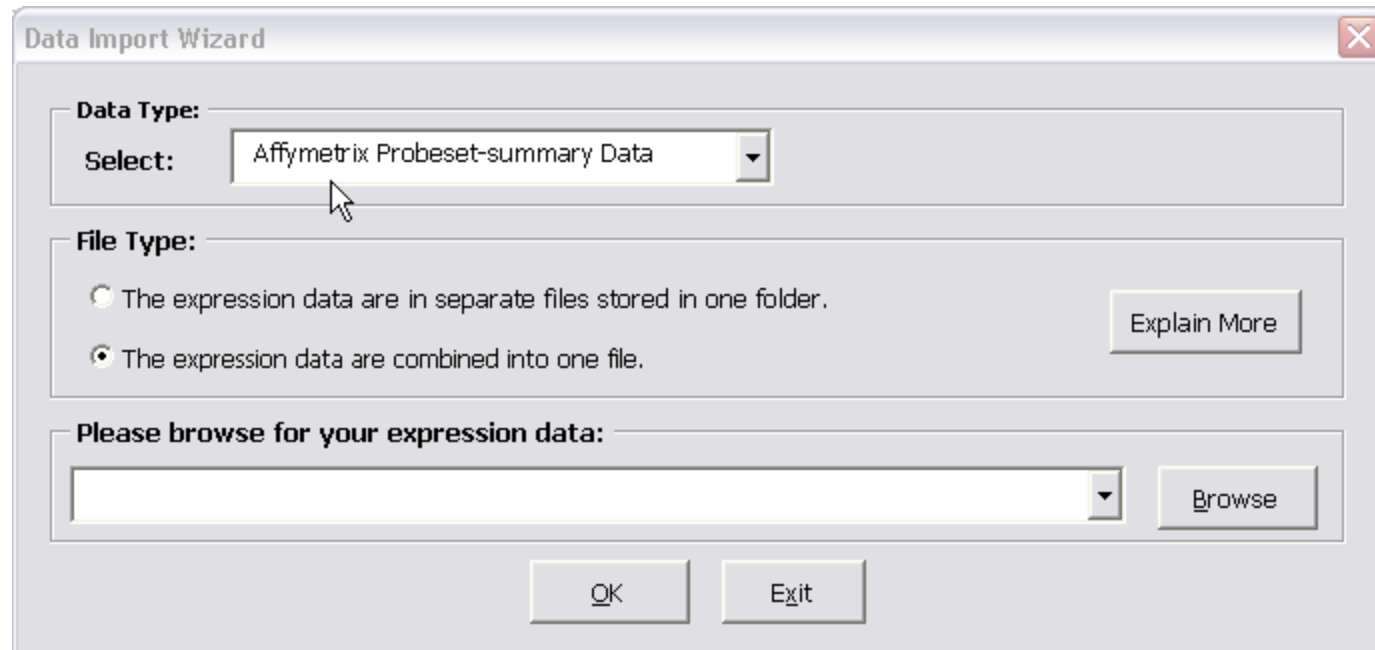


The screenshot shows the 'Data Import Wizard' dialog box. It has a title bar with a close button (X). The dialog is divided into three main sections. The first section, 'Data Type:', contains a 'Select:' label and a dropdown menu currently showing 'Affymetrix Probeset-summary Data'. A mouse cursor is pointing at the dropdown arrow. The second section, 'File Type:', contains two radio button options: 'The expression data are in separate files stored in one folder.' (unselected) and 'The expression data are combined into one file.' (selected). To the right of these options is an 'Explain More' button. The third section, 'Please browse for your expression data:', contains a text input field and a 'Browse' button. At the bottom of the dialog are 'OK' and 'Exit' buttons.

# [Hands-on instructions]

## [Importing Pomeroy Data set]


• **Browse** for the following file which is also in the **Pomeroy** folder inside the **BRB-ArrayTools Class** folder which is on the **Desktop**: **Dataset\_A2\_multiple\_tumor\_samples.txt** and then click **OK**.

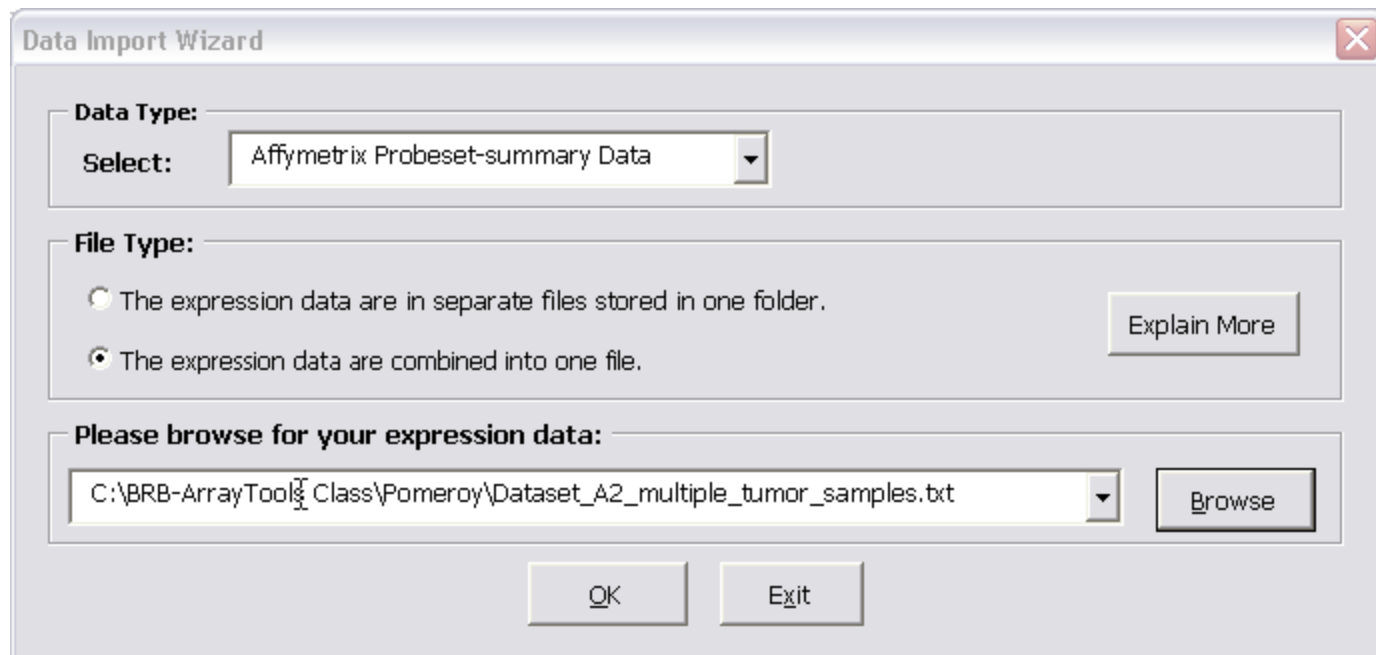


# [Hands-on instructions]

[Importing Pomeroy Data set]

- Click “yes” to the following question on number of arrays.

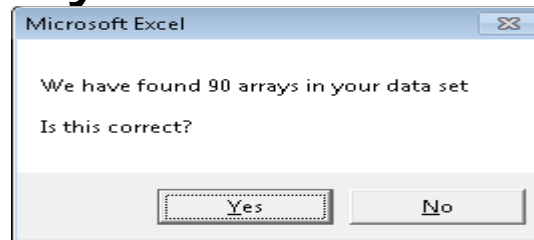


- Click “No” to the question about log transformation.

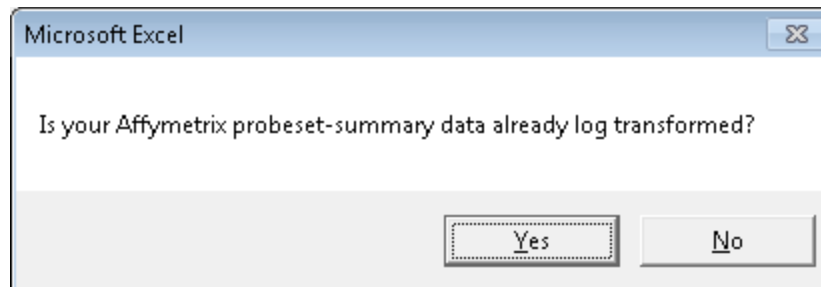

# [Hands-on instructions]

[Importing Pomeroy Data set]


- Select the chip type as “HuGeneFL Genome Array”

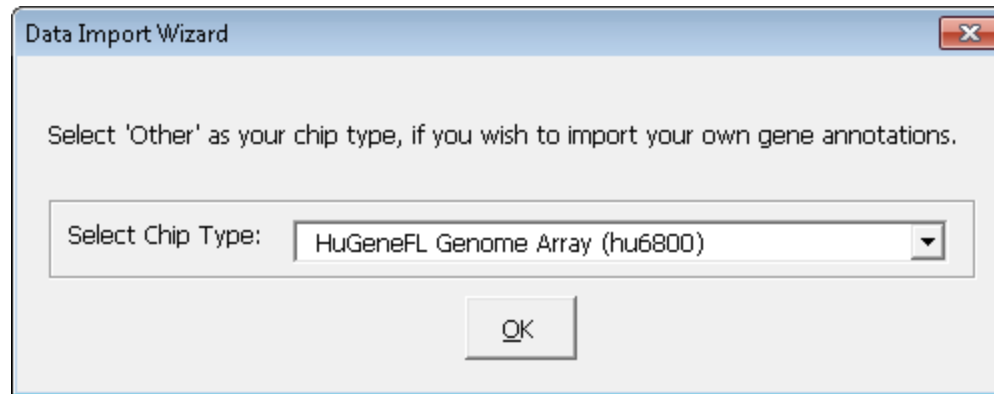

# [Hands-on instructions]

## [Importing Pomeroy Data set]


- **Browse** for the following file in the `Pomeroy` folder inside the `BRB-ArrayTools class` folder which is on the **Desktop**:  
“ExpDescrMedulo.xls” and click “Next”.

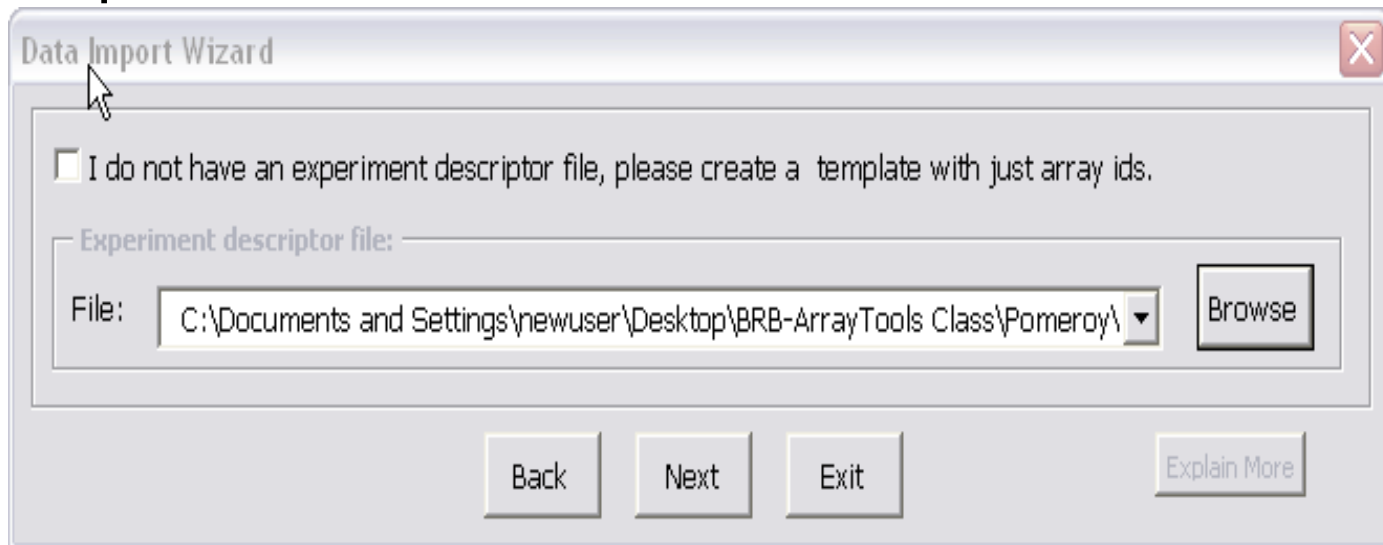

# [Hands-on instructions]

[Importing Pomeroy Data set]

- Save the Project in the folder “Pomeroy-Project”.

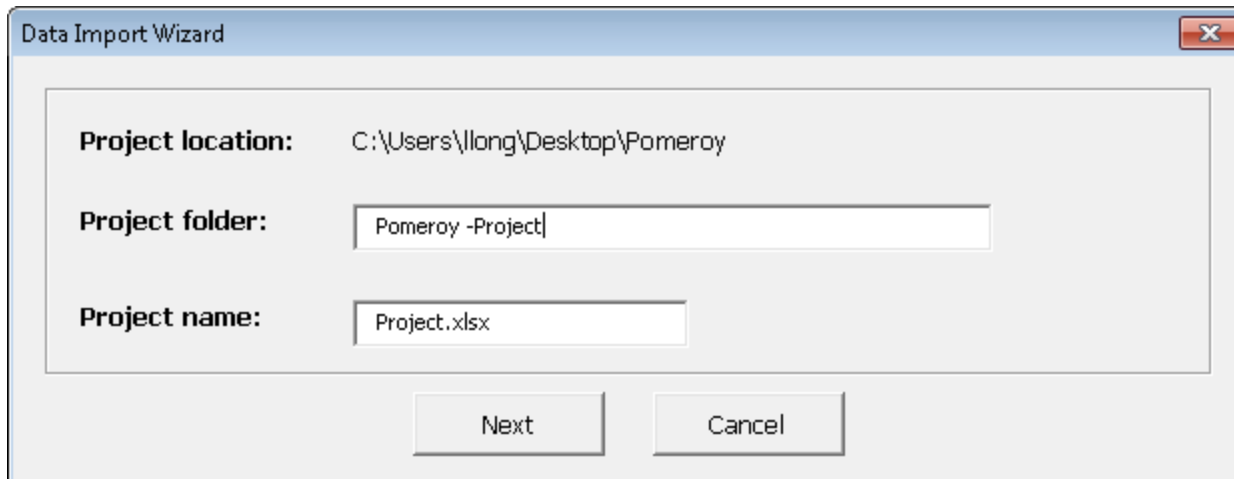

The screenshot shows a "Data Import Wizard" dialog box with a standard Windows window border. The title bar is light blue with a close button (X) on the right. The main area is light gray and contains three labeled fields: "Project location:" with the text "C:\Users\llong\Desktop\Pomeroy", "Project folder:" with a text box containing "Pomeroy -Project", and "Project name:" with a text box containing "Project.xlsx". At the bottom, there are two buttons: "Next" and "Cancel".


<b>Project location:</b>	C:\Users\llong\Desktop\Pomeroy
<b>Project folder:</b>	<input type="text" value="Pomeroy -Project"/>
<b>Project name:</b>	<input type="text" value="Project.xlsx"/>

Next Cancel

# [Hands-on instructions]

[Importing Pomeroy Data set]

- The progress bar will indicate that the project is collating.

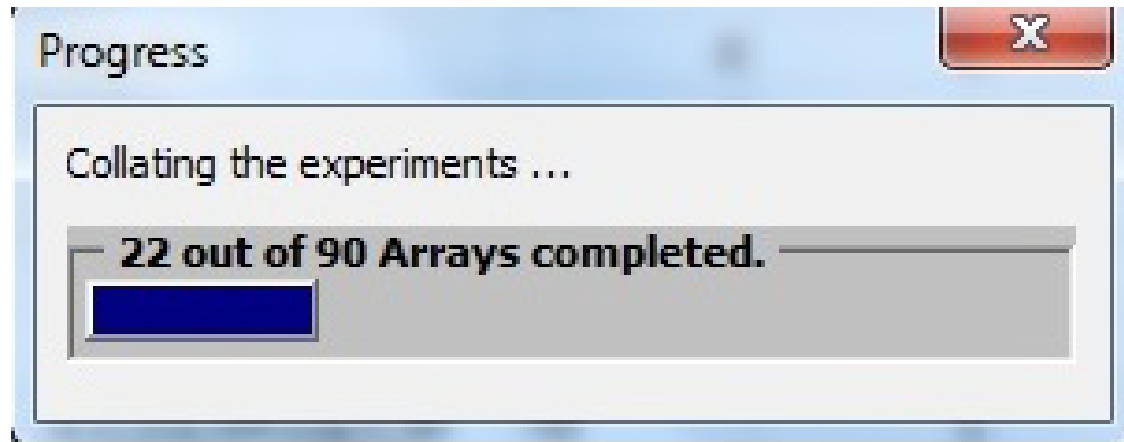



# Filtering and Normalization

Refilter, normalize and subset the data

1. Spot filters   2. Normalization   3. Gene filters

Background adjustment is performed before the intensity filtering and the averaging of replicate spots is done on filtered data.

☐ Apply background adjustment.

☒ **Intensity Filter:**

☐ EXCLUDE the spot if the intensity is below the minimum.

☒ THRESHOLD the intensity at the minimum value if the intensity is below the minimum.

Intensity minimum:

☐ **Detection Call**


EXCLUDE the probeset if the Detection Call contains:

Any one of the following values (comma-separated): A,M,P,No Cal


☐ Average the replicate spots within an array.

☐ Use a common reference design.

# When importing is completed

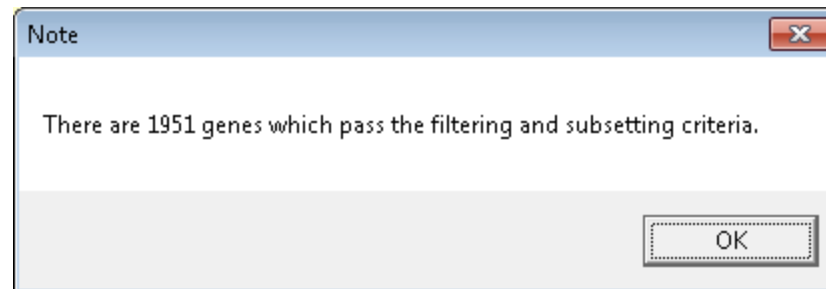

# Annotation

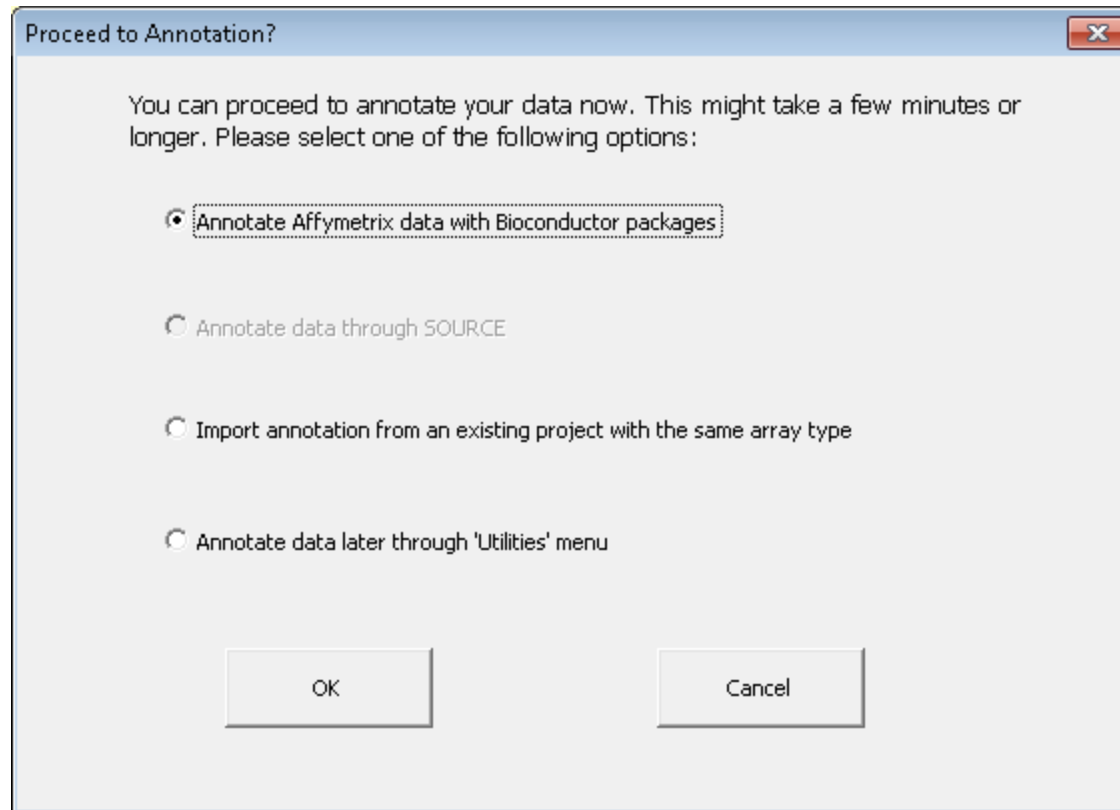


# Collated project workbook

## Overview

- The **collated project workbook** is the primary data object on which future analyses are run
- The collated project workbook is located inside the **project folder**, which by default is located inside the folder where the original input data is located.
- The project folder may also contain some other folders: **BinaryData**, **Annotations**, **Output**, and **Genelists**.
- The **BinaryData** and **Annotations** folders should NOT be altered by users. These are used for internal purposes.
- The **output** folder will contain the output of all subsequent analyses.
- A **Genelists** folder may also be created, and may contain genelists to be used for subset analyses.

# The collated project workbook

- This is the primary data object on which future analyses are run.
- Contains three primary worksheets:
  1. Experiment descriptors (may edit this to specify analyses)
  2. Gene identifiers
  3. Filtered log ratio (or Filtered log intensity)
- Additional results worksheets which may be automatically added:
  1. Gene annotations (obtained by running the menu item:  
**Utilities → Annotate data →**
  2. Cluster analysis results

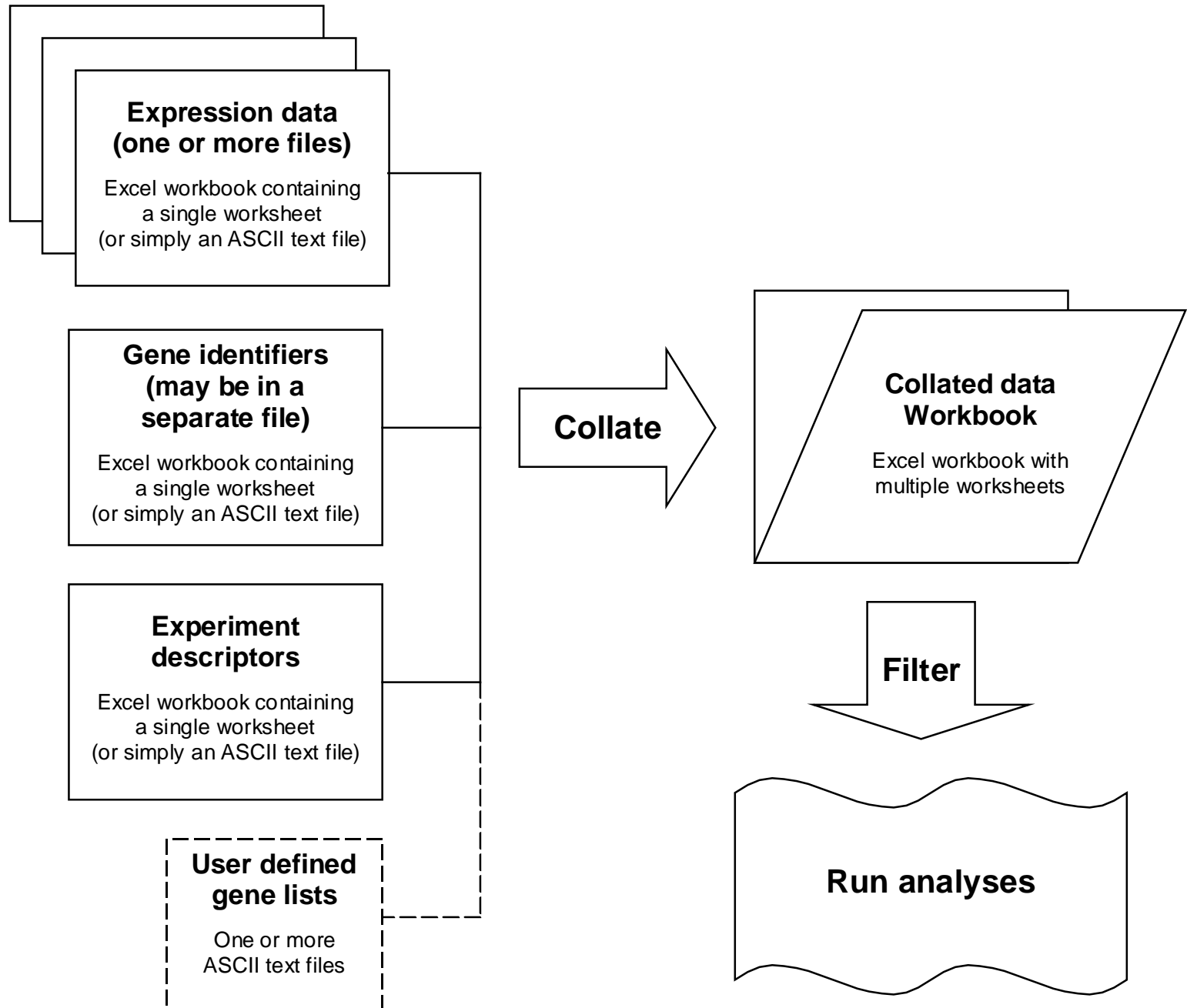

# The collated project workbook

## Experiment descriptor sheet

Create experiment descriptor variables which can be used to guide and specify the analyses.

Home Insert Page Layout Formulas Data Review View Add-Ins																
ArrayTools ▾																
CGHTools ▾																
Menu Commands																
A1 Array																
	A	B	C	D	E	F	G	H	I	J	K	L	M	N	O	P
1	Array	Dx	Medulo T	Medulo St	Sex	Age at Dx	Survival (r	SurvStatu	Chemo	AT RT Site	M Stage	SurvStatu	Age Code	Medulo vs	Glio vs	PNRhabdo
2	Brain_MD	Medullob	Classic	T4M1	M	8m	11 D		V,C,Cx,VP		>0	1	0	Medulloblastoma		
3	Brain_MD	Medullob	Classic	T2M0	M	8yr10m	5 D		V,C,Cx,VP		0	1	1	Medulloblastoma		
4	Brain_MD	Medullob	Classic	T3M0	M	6yr	7 D		V,C,Cx		0	1	1	Medulloblastoma		
5	Brain_MD	Medullob	Classic	T3M3	M	5yr 3m	7 D		V,C,Cx,VP		>0	1	1	Medulloblastoma		
6	Brain_MD	Medullob	Classic	M3	M	38yr 2m	7 D		V,C		>0	1	2	Medulloblastoma		
7	Brain_MD	Medullob	Classic	T4M0	F	7m	9 D		V,C,Cx		0	1	0	Medulloblastoma		
8	Brain_MD	Medullob	Classic	T1M0	M	6yr 5m	14 D		V,C,Cx		0	1	1	Medulloblastoma		
9	Brain_MD	Medullob	Classic	T3bM1	M	6yr 1m	16 D		V,C,Cx		>0	1	1	Medulloblastoma		
10	Brain_MD	Medullob	Classic	M0	M	8yr	18 D		V,C,Cx,VP		0	1	1	Medulloblastoma		
11	Brain_MD	Medullob	Classic	M0	M	3yr 10m	18 D		V,C,Cx		0	1	0	Medulloblastoma		
12	Brain_MD	Medullob	Classic	T2M1	M	8yr 2m	19 D		V,C,Cx,VP,Ca,T,M		>0	1	1	Medulloblastoma		
13	Brain_MD	Medullob	Classic	M0	F	3yr 9m	25 D		V,C,Cx		0	1	0	Medulloblastoma		
14	Brain_MD	Medullob	Classic	T3M3	M	14yr 5m	26 D		V,C,Cx		>0	1	1	Medulloblastoma		
15	Brain_MD	Medullob	Desmopla	M0	M	6yr 3m	33 D		V,C,CC		0	1	1	Medulloblastoma		
16	Brain_MD	Medullob	Desmopla	T2MO	F	11yr 7m	38 D		V,C,Cx,VP		0	1	1	Medulloblastoma		
17	Brain_MD	Medullob	Desmopla	T3M3	F	11yr 5m	39 D		V,C,VP		>0	1	0	Medulloblastoma		
18	Brain_MD	Medullob	Classic	T3bM3	F	3yr 3m	39 D		V,C,Cx		>0	1	0	Medulloblastoma		
19	Brain_MD	Medullob	Classic	T2M3	M	4yr 4m	42 D		V,C,Cx		>0	1	0	Medulloblastoma		
20	Brain_MD	Medullob	Classic	M2	F	26yr 1m	65 D		V,C,Cx,VP		>0	1	2	Medulloblastoma		
21	Brain_MD	Medullob	Classic	T3bM0	M	20yr 6m	92 D		V,C		0	1	2	Medulloblastoma		
22	Brain_MD	Medullob	Classic	T2M0	F	23yr 3m	102 D		V,C		0	1	2	Medulloblastoma		
23	Brain_MD	Medullob	Desmopla	M0	F	5yr 7m	24 A		V,C,CC		0	0	1	Medulloblastoma		
24	Brain_MD	Medullob	Desmopla	T4M0	M	1yr 4m	25 A		V,C,Cx		0	0	0	Medulloblastoma		

# The collated project workbook

## Gene identifier sheet

Contains gene identifiers provided by the user during collation.

A1 ProbeSet						
	A	B	C	D	E	F
1	ProbeSet	Name	Symbol	Entrez	Defined_Genelists	Filter
20	hum_alu_at					TRUE
39	AFFX-HUN	signal tran	STAT1	6772	Apoptotic Signaling in Respo	TRUE
43	AFFX-HUN	glyceralde	GAPDH	2597	Downregulated of MTA-3 in E	TRUE
44	AFFX-HUN	glyceralde	GAPDH	2597	Downregulated of MTA-3 in E	TRUE
45	AFFX-HUN	glyceralde	GAPDH	2597	Downregulated of MTA-3 in E	TRUE
46	AFFX-HSA	actin, bet	ACTB	60	Chromatin Remodeling by hS	TRUE
47	AFFX-HSA	actin, bet	ACTB	60	Chromatin Remodeling by hS	TRUE
48	AFFX-HSA	actin, bet	ACTB	60	Chromatin Remodeling by hS	TRUE
52	AFFX-M27830_5_at					TRUE
55	AFFX-HSA	actin, bet	ACTB	60	Chromatin Remodeling by hS	TRUE
58	AFFX-HUN	glyceralde	GAPDH	2597	Downregulated of MTA-3 in E	TRUE
64	AB000220	sema dom	SEMA3C	10512	Axon guidance	TRUE
68	AB000460	family wit	FAM193A	8603		TRUE
70	AB000464	NOP14 an	NOP14-AS	317648		TRUE
73	AB000468	ring finger	RNF4	6047		TRUE
79	AB001106	glia matur	GMFB	2764		TRUE
80	AB001325	aquaporin	AQP3	360		TRUE
82	AB002315	KIAA0317	KIAA0317	9870		TRUE
86	AB002380	Rho guani	ARHGEF12	23365	Axon guidance, Regulation o	TRUE
87	AB002382_at					TRUE
90	AB003102	proteasom	PSMD11	5717	Proteasome	TRUE
91	AB003103	proteasom	PSMD12	5718	Proteasome	TRUE
92	AB003177	proteasom	PSMD9	5715		TRUE
93	AB003698	cell divis	CDC7	8317	Cell cycle	TRUE
94	AB004884	tousled-li	TLK2	11011		TRUE
99	AC000064	GATA zinc	GATAD1	57798		TRUE



# The collated project workbook

Filtered log ratio or log intensity sheet

View the matrix of log-expression data with data filters applied.

3.5 / 34033584594 /														
	A	B	C	D	E	F	G	H	I	J	K	L	M	N
	Display the data	Brain_M D_1	Brain_M D_2	Brain_M D_3	Brain_M D_4	Brain_M D_5	Missing	P-Value	Rank	Variance	Num 1.5- Fold	Absent (Affy)	Filter	
1	ProbeSet													
20	hum_alu_at	14.27817	14.44255	15.29377	14.9108	16.8653	0				57	3	TRUE	
39	AFFX-HUMISGF3A/M97	8.927432	8.714246	8.096399	7.865125	10.75989	0				53	22	TRUE	
43	AFFX-HUMGAPDH/M33	14.27731	10.03617	14.0167	13.4024	13.36646	0				18	0	TRUE	
44	AFFX-HUMGAPDH/M33	14.1479	10.66977	14.57995	13.45807	13.64926	0				26	1	TRUE	
45	AFFX-HUMGAPDH/M33	14.28507	13.62297	15.00179	14.45779	14.12614	0				28	0	TRUE	
46	AFFX-HSAC07/X00351_5	14.13774	9.240791	14.31056	13.54756	14.78994	0				34	2	TRUE	
47	AFFX-HSAC07/X00351_1	14.61301	11.43411	14.84308	14.13051	15.13767	0				36	2	TRUE	
48	AFFX-HSAC07/X00351_3	14.08011	14.59584	14.15988	13.69802	14.62119	0				26	0	TRUE	
52	AFFX-M27830_5_at	10.9794	9.914385	11.43107	8.792543	11.89633	0				78	38	TRUE	
55	AFFX-HSAC07/X00351_3	8.794507	10.46557	10.64404	9.315951	11.88684	0				54	5	TRUE	
58	AFFX-HUMGAPDH/M33	7.084365	9.240791	8.576831	8.371586	9.057992	0				66	44	TRUE	
64	AB000220_at	9.557537	7.882643	6.681361	7.042521	9.971543	0				61	34	TRUE	
68	AB000460_at	10.15234	10.83368	10.6047	11.35548	10.23002	0				32	1	TRUE	
70	AB000464_at	8.536878	10.04576	8.547344	9.687808	10.29692	0				53	27	TRUE	
73	AB000468_at	9.438828	7.761551	8.468632	3.445586	9.108524	0				52	15	TRUE	
79	AB001106_at	7.702686	9.647458	10.19205	6.994023	8.816983	0				55	14	TRUE	
80	AB001325_at	8.140219	10.37069	9.846106	7.858368	10.03342	0				37	12	TRUE	
82	AB002315_at	6.790634	7.066089	6.842642	3.445586	7.569856	0				54	36	TRUE	
86	AB002380_at	8.591326	8.169925	6.870123	6.662817	9.562243	0				68	23	TRUE	
87	AB002382_at	10.56322	9.569856	8.764255	8.778295	9.473705	0				39	32	TRUE	
90	AB003102_at	10.52527	9.098032	9.39959	7.803138	10.15861	0				54	7	TRUE	
91	AB003103_at	6.743329	7.721099	8.689096	6.994023	8.523562	0				54	18	TRUE	
92	AB003177_at	8.568888	8.430452	7.814628	7.32333	9.092757	0				45	26	TRUE	
93	AB003698_at	6.821331	8.144658	7.778827	5.515975	8.754888	0				51	39	TRUE	

# The collated project workbook

## Gene annotations worksheet (Optional)

Contains gene annotations which were automatically downloaded from the Affymetrix or SOURCE database using the annotations tool.

	A	B	C	D	E	F	G	H	I	J	K
1	ProbeSet (Double-click)	Name	Accessi	UGClus	Symbol	Entrezl	Chrom	Cytoba	GO	Filter	
20	<a href="#">hum_alu_at</a>		U14573							TRUE	
39	<a href="#">AFFX-HUN</a>	signal tran	M97935	Hs.642990	STAT1	6772	2	2q32.2	Biological	TRUE	
43	<a href="#">AFFX-HUN</a>	glyceralde	M33197	Hs.544577	GAPDH	2597	12	12p13	Biological	TRUE	
44	<a href="#">AFFX-HUN</a>	glyceralde	M33197	Hs.544577	GAPDH	2597	12	12p13	Biological	TRUE	
45	<a href="#">AFFX-HUN</a>	glyceralde	M33197	Hs.544577	GAPDH	2597	12	12p13	Biological	TRUE	
46	<a href="#">AFFX-HSA</a>	actin, beta	X00351	Hs.520640	ACTB	60	7	7p22	Molecular	TRUE	
47	<a href="#">AFFX-HSA</a>	actin, beta	X00351	Hs.520640	ACTB	60	7	7p22	Molecular	TRUE	
48	<a href="#">AFFX-HSA</a>	actin, beta	X00351	Hs.520640	ACTB	60	7	7p22	Molecular	TRUE	
52	<a href="#">AFFX-M27830_5_at</a>		M27830							TRUE	
55	<a href="#">AFFX-HSA</a>	actin, beta	X00351	Hs.520640	ACTB	60	7	7p22	Molecular	TRUE	
58	<a href="#">AFFX-HUN</a>	glyceralde	M33197	Hs.544577	GAPDH	2597	12	12p13	Biological	TRUE	
64	<a href="#">AB000220</a>	sema dom	AB000220	Hs.269109	SEMA3C	10512	7	7q21-q31	Biological	TRUE	
68	<a href="#">AB000460</a>	family wit	AB000460	Hs.652364	FAM193A	8603	4	4p16.3	Cellular C	TRUE	
70	<a href="#">AB000464</a>	NOP14 an	AB000464	Hs.398178	NOP14-AS	317648	4	4p16.3		TRUE	
73	<a href="#">AB000468</a>	ring finger	AB000468	Hs.740360	RNF4	6047	4	4p16.3	Molecular	TRUE	
79	<a href="#">AB001106</a>	glia matur	AB001106	Hs.151413	GMFB	2764	14	14q22.2	Molecular	TRUE	
80	<a href="#">AB001325</a>	aquaporin	AB001325	Hs.234642	AQP3	360	9	9p13	Biological	TRUE	
82	<a href="#">AB002315</a>	KIAA0317	AB002315	Hs.730659	KIAA0317	9870	14	14q24.3	Molecular	TRUE	
86	<a href="#">AB002380</a>	Rho guani	AB002380	Hs.24598	ARHGEF12	23365	11	11q23.3	Molecular	TRUE	
87	<a href="#">AB002382_at</a>		AB002382						Biological	TRUE	
90	<a href="#">AB003102</a>	proteasom	AB003102	Hs.595584	PSMD11	5717	17	17q11.2	Biological	TRUE	
91	<a href="#">AB003103</a>	proteasom	AB003103	Hs.592689	PSMD12	5718	17	17q24.2	Biological	TRUE	
92	<a href="#">AB003177</a>	proteasom	AB003177	Hs.131151	PSMD9	5715	12	12q24.31-	Biological	TRUE	
93	<a href="#">AB003698</a>	cell divis	AB003698	Hs.533573	CDC7	8317	1	1p22	Biological	TRUE	
94	<a href="#">AB004884</a>	tousled-li	AB004884	Hs.445078	TLK2	11011	17	17q23	Biological	TRUE	

## Part IV:

# Data filtering and normalization options

# [Hands-on instructions]

[Data filtering-Pomeroy]

1. Click on **ArrayTools** → **Re-Filter, normalize and subset the data.**
2. Click on the four buttons **Spot filter, Normalization, Gene filter and Gene Subset** at the TOP of the form, to see the available options and view the current settings applied on the dataset.
3. By clicking “OK” the default filtering and normalization is performed on the data set.

# Data filtering options

## Spot filter

- Intensity filter: May filter out spots with low intensity or threshold low intensity.
- Spot size filter: May filter out spots whose sizes less than a specified value.
- Spot flag filter: May filter out spots outside a range or containing specific values.
- Detection Call: Exclude a probeset if the Detection call value is “A”, “M”, “P” or “No Call”-specific for Affymetrix data
- Option of averaging replicate spots

# Data filtering options

## Normalization and truncation

- Normalization and truncation steps are applied *after* data has been spot-filtered, but *before* screening out genes
- Arrays are normalized before outlying expression levels are truncated.
- Purpose of truncation is primarily to prevent extremely large ratios from being formed by small denominators in dual-channel data. The truncation option is useful if the dual-channel intensities have not been thresholded.

# Data filtering options

## Data transformation options

- Normalization:

For single-channel data: Default option is quantile normalization.

For dual-channel data: Default option is median normalization.

- Truncation: Truncate extreme values (large log-intensities for single-channel data, or large absolute log-ratios for dual channel data)

# Data filtering options

## Gene filters: Gene variation

- Fold-change filter: Specify a minimum percentage of log-expression values which must meet a specified fold-change criteria
- Log-ratio (or log-intensity) variation filter:  
Screen genes which do not vary much over the set of samples:
  1. Significance criterion compares the variance of each gene against the “average” gene
  2. Percentile criterion screens a specified percentage of genes with smallest variance



# Data filtering options

## Gene filters: Gene quality

- Missing value filter: Screens out genes which contain too many missing values over the set of samples
- Percent absent filter: For Affymetrix data, can filter out a probeset if too many expression values had an Absent call
- Minimum Intensity: This option is only available for single channel data. It filters out genes whose specified percentile normalized log intensity is less than the log of the user defined value.

# Data filtering options

## Gene subsets

- Select genelists for analysis: User may subset the data by selecting one or more genelists to INCLUDE or EXCLUDE. If more than one genelists is selected, then the UNION of all genes on those genelists will be used.
- Specify gene labels to exclude: User may exclude genes based on gene identifier labels. For example, all genes with “Empty” in the gene description field may be excluded.
- CAUTION: Gene subsetting is applied globally to the entire dataset, not just to a specific analysis.
- Probe reduction: Reduce multiple probe sets per gene by choosing the most variably expressed or the maximally expressed probe/probeset.

# Spot Filtering options- Dual Channel data

Refilter, normalize and subset the data

1. Spot filters   2. Normalization   3. Gene filters

Background adjustment is performed before the intensity filtering and the averaging of replicate spots is done on filtered data.

☐ Apply background adjustment.

☒ **Intensity Filter:**

- ☐ EXCLUDE the spot if BOTH intensities are below the minimum.
- ☐ EXCLUDE the spot if AT LEAST ONE of the two intensities is below the minimum.
- ☒ EXCLUDE the spot if BOTH intensities are below the minimum. If only ONE intensity is below the minimum, increase it to the minimum.

Red minimum:

Green minimum:

☐ **Spot Flag Filter:**

EXCLUDE the spot if the Spot Flag contains:

☒ Numeric values outside the range:  to

☐ Any of the following values:

List of values, separated by commas:

☐ **Spot Size Filter:**

EXCLUDE the spot if the Spot Size is less than:

☐ Average the replicate spots within an array.

☐ Use a common reference design.

# Normalization options- Dual channel data

Refilter, normalize and subset the data

1. Spot filters   2. Normalization   3. Gene filters   4. Gene subsets

A log base 2 transformation is applied to the data before the arrays are normalized.

☒ **Normalize (center) each array:**

☒ Using median over entire array   ☐ Using median with print-tip group  
☐ Using lowess smoother   ☐ Using Lowess with print-tip group  
☐ Using median over housekeeping genes:

☐ HG-U133 (Affy)   ☐ HG-U95 (Affy)   ☐ HG-Focus (Affy)  
☒ Specify a set of housekeeping genes:

☐ **Truncate large intensity ratios (and inverse ratios):**

Truncate intensity ratios greater than:

# Part V:

## Scatterplot tools

# Scatterplot tools

- Scatterplot of array vs. array: Plots intensity, geometric mean of the red and green intensities, and intensity ratio on log-scale.
- Scatterplot of phenotype averages: Plots averages over experiment classes
- Online demo  
<http://linus.nci.nih.gov/PowerPointSlides/Scatterplot.wmv>

# [Optional: Hands-on instructions]

## [Scatterplot of phenotype averages]

1. Now click on **ArrayTools** → **Graphics** → **Scatterplot** → **Phenotype averages**.
2. Select the variable **Dx** as the phenotype class to average over, and then click **OK**.
3. This launches a 2-D and 3-D scatter plot.
4. Right click on the 2-D plot to modify scatter plot properties, select up/down regulated genes as well as link genes in other plots.

# [Optional: Hands-on instructions]

[Scatterplot of experiment v. experiment-Pomeroy Data]

1. Click on **ArrayTools** → **Graphics** → **Scatterplot** → **Array vs. Array**.
2. Select **Log(Intensity)** for the **Brain\_MD\_1** experiment for the X-values and **Log(Intensity)** for the **Brain\_MD\_MGlio\_1** experiment as Y-values.
3. Select “2” as the number of panels.
4. Click “OK”. Then, right click on the plot to change scatterplot properties, select up/down regulated genes etc.



# Part VI:

## Class Comparison

# Class Comparison

## Experimental design

Step1

Step2

Experimental design:

Column defining classes:

☒ Unpaired samples:

☐ Block by:

☐ Average over replicates of:

☐ Paired samples:

Pair samples by:

# Class comparison tool

1. Enter the class column from the 'Experiment descriptor' worksheet that defines the classes for the samples.
2. Specify if this is a paired or un-paired analysis. An analysis is said to be paired if for example, you have the same sample from a patient before and after a treatment. You then need a column in the experiment descriptor worksheet that will contain identical values for pair of arrays.
3. If this is an unpaired analysis, do you have a [blocking factor?](#)
4. If this is an unpaired analysis, do you have an [replicates](#) you want to average across?

# Class comparison tool

## Blocking Factor

Experimental designs containing a blocking factor can be performed by specifying which column in the Experiment descriptor worksheet contains a blocking variable. When selected, the influence of the blocking variable is taken into consideration when analyzing the differences between classes.

Examples of variables that may be considered as Blocking factors:

- Clinical Site for patient data
- Print set for cDNA spotted arrays
- Batch of arrays


# Average over replicates

- If multiple arrays have been performed using the same sample RNA then an average of these replicates should be used instead of the individual arrays in the analysis.
- In the 'experiment descriptor' worksheet, there should be column containing sample ids for these arrays.
- Arrays that contain the identical values of the sample id variable are considered as replicates and will be averaged in the analysis.

# Class comparison tool

## Random variance option

- The random variance test has more power because the “average” variance in the denominator adds degrees of freedom for the test statistic.
- Should be used for small sample sizes.
- Dialog option:

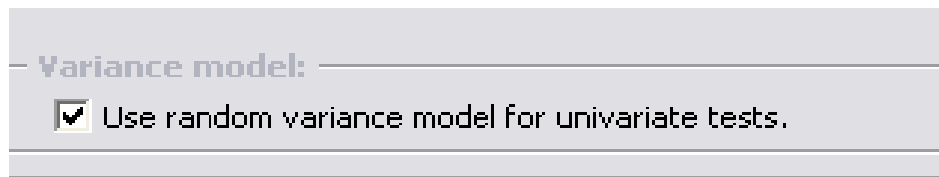

# Find genes lists determined by:

**Find gene lists determined by:**

☒ Significance threshold of univariate tests: 0.001

☐ Multivariate permutation tests

Maximum proportion of false discoveries: 0.1

Confidence level (between 0 and 100%): 80

☐ Threshold of the local false discovery rate:  
0.05

# Class comparison tool

## Univariate significance test

- Compute the univariate p-value for each gene, and sort list of genes by smallest p-value.
- In the univariate setting (i.e., testing significance of one gene at a time), the p-value is defined to be the probability of obtaining a false positive result.
- However, once a list of univariately significant genes is found, it is not clear how many of those genes are false positives.



# [Hands-on instructions]

## [Class comparison – univariate significance threshold]

1. Using the Pomeroy data, run the Class Comparison tool by clicking on **ArrayTools → Class comparison → Between groups of arrays**.
2. Select the **Medulo vs Glio** variable as the column defining the classes. Select the **Random variance model** option, and select the **Significance threshold of univariate tests: 0.001**.
3. Leave all other options at default levels. Now click **OK** on the main dialog to launch the analysis.
4. You will see a DOS window appear in your Windows Task Bar at the bottom of your screen. If you click on the DOS window, you can monitor the analysis running inside the DOS window.
5. When the analysis has completed, it will automatically open up an HTML file which displays the output.

# Class comparison tool

## Multivariate permutation test

Find gene lists determined by:

☐ Significance threshold of univariate tests: 0.001

☒ Multivariate permutation tests

Maximum proportion of false discoveries: 0.1

Confidence level (between 0 and 100%): 80

☐ Threshold of the local false discovery rate: 0.05

# Class comparison tool

## Multivariate permutation test

- In the multivariate setting (i.e., when testing many genes for significance at the same time), ask the question: What p-value cutoff should I use to guarantee that 90% of the time, I get less than  $P$  proportion of false positives (where  $P$  is specified by the user)?
- To answer this question, we compute the permutation distribution of the p-value cutoffs for which we would get  $P$  proportion of false positives.
- The output tells us how far down the list we would be able to go in order to be assured (with a certain confidence) of getting less than  $P$  proportion of false positives.

# [Hands-on instructions]

[Class comparison – Restricting proportion of false positives]

1. Using the Pomeroy data, run the Class Comparison tool by clicking on **ArrayTools** → **Class comparison** → **Between groups of arrays**.
2. Select the **Medulo vs Glio** variable as the column defining the classes. Select the **Random variance model** option, and select the **Restriction on proportion of false discoveries** with **maximum proportion** = 0.1 and 90% **Confidence level**.
3. Click on the **options** and change the name of the **output** folder to “ClassComparisonMPT”
4. Leave all other options at default levels. Now click **OK** on the main dialog to launch the analysis.
5. When the analysis has completed, it will automatically open up an HTML file which displays the output.

# Class comparison tool

## Local False Discovery Rate (LFDR)

Find gene lists determined by:

☐ Significance threshold of univariate tests:

0.001

☐ Multivariate permutation tests

Maximum proportion of false discoveries:

0.1

Confidence level (between 0 and 100%):

80

☒ Threshold of the local false discovery rate:

0.05

# Class comparison tool

## Local False Discovery Rate

- An empirical Bayes approach for two-class or paired samples microarray data.
- Find differentially expressed genes by controlling for local false discovery rate.

*Efron, B., Tibshirani, R., Storey, J.D., and Tusher, V.  
Empirical Bayes analysis of a microarray experiment. Journal  
of the American Statistical Association, 96(456):1151-1160,  
2001.*

# [Hands-on instructions]

[Class comparison – Threshold LFDR]

1. Using the Pomeroy data, run the Class Comparison tool by clicking on **ArrayTools** → **Class comparison** → **Between groups of arrays**.
2. Select the **Medulo vs Glio** variable as the column defining the classes. Select the **Random variance model** option, and select the **Threshold of Local Discovery Rate = 0.05**.
3. Click on the **options** and change the name of the **output** folder to “ClassComparisonLFDR”
4. Leave all other options at default levels. Now click **OK** on the main dialog to launch the analysis.
5. When the analysis has completed, it will automatically open up an HTML file which displays the output.

# Class comparison

## Significance Analysis of Microarrays (SAM)

- SAM is another popular method for false discovery control, which controls the **average** proportion of false discoveries rather than the **probability** of a given number or proportion of false discoveries.
- It is a slightly less stringent control than the multivariate permutation test for controlling false discoveries used in the other class comparison tools, but is included in BRB-ArrayTools because of its popularity.



# [Hands-on instructions]

[Significance Analysis of Microarrays – Pomeroy data]

1. Still using the Pomeroy data, run the SAM tool by clicking on **ArrayTools → Class comparison → Significance Analysis of Microarrays (SAM)**.
2. Again, select the **Medulo vs Glio** variable as the column defining the classes, select the **90<sup>th</sup> percentile** option, and leave all other parameters at default levels.
3. Check the option to perform **Gene ontology Observed vs Expected analysis**.
4. Now click **OK** to exit the options dialog, and click **OK** on the main dialog to launch the analysis.

# Part VII:

## Gene Set Expression Comparison

# Gene set Expression Comparison

- Allows users to find significant **sets** of genes rather than just significant genes.
- For the **Gene Ontology comparison**, all Gene Ontology classes that are represented in the data are tested for significance.
- For **Pathway Comparison**, all the pathways that are represented in the data are tested. For Human, the BioCarta or KEGG pathways are tested and for mouse, the BioCarta pathways are compared. Additionally, Broad/MIT pathways can be downloaded to be used in analyses.
- For the **User Gene Lists comparison**, the user can select specific genelists that the user would like to test for significance.
- Transcription factor target gene lists and microRNA target genelists have been added to the Gene List comparison tool.

# Gene Set Expression Comparison

- Compute p-value of differential expression for each gene in the gene set( $k$ =number of genes)
- Compute a summary ( $S$ ) of these p-values
- Determine whether the summary test ( $S$ ) is more extreme than would be expected from a random sample of “ $k$ ” genes on that platform.
- Two types of summaries provided:
  - Average of log p-values
  - Kolmogorov-Smirnov statistic.

# Efron-Tibshirani's GSA maxmean test

- Tests the null hypothesis that for a gene set the average degree of differential expression is greater than expected from a random set of genes.
- Uses the maxmean statistic as follows:
- Take the  $d_i$  scores for all the genes within a geneset.
- Set negative scores to 0 and compute 'avpos' as the average of the positive scores and zeros.
- Similarly set the positive scores to 0 and compute the 'avneg' as the averages of the negative scores and zeros.
- A gene set is scored 'avpos' if  $|avpos| > |avneg|$  or else the gene set is scored 'avneg'

# [Hands-on instructions]

[Class Comparison – Pathway Comparison: Pomeroy data]

1. On the Pomeroy data, run the Class Comparison tool by clicking on **ArrayTools → Class comparison → Gene set Expression Comparison**.
2. Select the **Medulo vs Glio** variable as the column defining the classes. Select the **Random variance model** option and **Pathways**, and leave all other options at default levels. Now click **OK** on the main dialog to launch the analysis.
3. You will see a DOS window appear in your Windows Task Bar at the bottom of your screen. If you click on the DOS window, you can monitor the analysis running inside the DOS window.
4. When the analysis has completed, it will automatically open up an HTML file which displays the output.

# Quantitative trait tool

- Selects genes which are univariately correlated with a quantitative trait such as age or time point.
- Controls number and proportion of false discoveries in entire list: uses a multivariate permutation test which takes advantage of the correlation among genes.
- Produces a gene list which can be used for further analysis.
- Produces chromosomal distribution and GO analysis if genes have already been annotated using the SOURCE database.

# Survival analysis tools

- Find Genes Correlated with Survival tool, selects genes which are univariately correlated with survival
- Controls number and proportion of false discoveries in entire list: uses a multivariate permutation test which takes advantage of the correlation among genes
- Produces a gene list which can be used for further analysis.
- Produces chromosomal distribution and GO analysis if genes have already been annotated using the SOURCE database.



# Survival Gene Set analysis

- This analysis tool finds sets of genes for which the expression levels are correlated to survival. Similar to the Gene Set Expression comparison tool, this tool can be used to analyze Gene Ontology categories, Pathways, micro RNA targets, transcription factor targets and user defined gene lists.
- The permutation p-values from the LS and KS statistics are computed.
- The HTML output lists the sets of genes and the associated p-values.

# **Part I: Clustering**

# Hierarchical clustering tools

- Clustering of genes and samples produces visual image plot of log-expression data, where ordering is determined by ordering of dendrogram
- Can compute measures to assess cluster reproducibility when clustering samples alone
- May cluster based on gene subsets rather than on the entire gene set
- Interface to Cluster 3.0 and TreeView originally produced by the Stanford group is also included, and allows for easy exportation of results.

# [Hands-on instructions]

[Cluster analysis – Pomeroy data]

1. Using the Pomeroy data set.
2. Run the cluster analysis by clicking on **ArrayTools** → **Clustering** → **Dynamic Heatmap Viewer**.
3. Click on the **Select gene subsets** button, and under **Select genes for analysis**, choose the **ClassComparison** genelist, and click **OK**.
4. Click “OK” on the main dialog page.

# [Hands-on instructions – cont'd]

[Cluster analysis – Pomeroy data]

10. Still with the Pomeroy data in front of you, click on the **ArrayTools** → **Clustering** → **Sample alone** menu item.
11. Select the **Compute the cluster reproducibility** option
12. Now click on the **Options** button, and choose **Dx** as the variable under **Label the experiments**.
13. Click **OK** to exit the options dialog, and click **OK** on the main dialog to launch the analysis.

# [Hands-on instructions – cont'd]

## [Cluster analysis – Pomeroy data]

- 14: The analysis will create a dendrogram plot of the hierarchical clustering of samples inside the **Cluster viewer** worksheet. You may then click the **Cut tree(# of cluster 3)** button to “cut the tree”, thereby defining clusters of samples from the dendrogram. After you have defined clusters of samples by “cutting the tree”, the analysis will be run in a DOS window which appears in your Windows Task Bar, and an HTML file containing the output will open up automatically once the computation is completed

# Cluster reproducibility

- Add perturbation noise to original data
- Re-cluster perturbed data to assess stability of original clusters
- Overall and cluster-specific measures
- Robustness (R) index measures the proportion of pairs of specimens within a cluster for which the members of the pair remain together in the re-clustered perturbed data
- Discrepancy (D) index measures the number of discrepancies (additions or omissions) comparing an original cluster to a best-matching cluster in the re-clustered perturbed data.

# Part II: MDS



# Multidimensional scaling

- Rotating scatterplot: Gives three-dimensional visualization of relationships between samples
- Global test of clustering in samples: Compares spatial distribution of data to white noise. Large deviation from Gaussian normal distribution indicates presence of clustering.

# [Hands-on instructions]

[Multidimensional scaling –Pomeroy data]

1. Still using the **Pomeroy** dataset, run the multidimensional scaling by clicking on **ArrayTools → Graphics -> 3D-Visualization of samples**.
2. Now choose **Dx** as the variable to **Color the rotating scatterplot**. click **OK** on the main dialog to launch the analysis.
3. An RGL device window will be launched, containing a scatterplot which can be rotated using arrow control buttons. Each point represents a sample, and points can be identified by brushing over them with your mouse.
4. A PowerPoint slide is automatically created, so that you can also launch the rotating scatterplot at a later point from PowerPoint.

# Part III:

## Class Prediction

# Analysis Wizard- Prediction

- Class Prediction
- PAM
- Top scoring pair plug-in
- Random Forest plug-in
- Binary Tree Prediction

# Components of Class Prediction

- C1. Feature(gene) selection
  - which genes will be included in the model.
- C2. Select model type.
  - choose prediction method (DLDA,CCP etc)
  - Fit the parameters for the model.
- C3. Evaluating the Classifier
  - Cross-validation

# C1. Gene Selection Criteria

- Selection of genes may be based on univariate significance criterion or univariate misclassification rate, and minimum fold-ratio of geometric means. The univariate misclassification rate criterion is available when there are only two classes. The option to optimize over a grid of alpha values.
- In addition, we have added the option to select genes using “gene pairs” by the “greedy pair” method –Bo & Jonassen

# Gene Selection Criteria

Gene selection

☒ Individual genes:

☒ Significant univariately at alpha level:

☐ Optimize over the grid of alpha-levels (and cross-validate optimization)

☐ With univariate misclassification rate below:

☐ With fold-ratio of geometric means between two classes exceeding:

☐ Gene pairs

Number of pairs selected by the "Greedy pairs" method:


☐ Recursive feature elimination

Number of features to be selected:

# C2. Class prediction Methods

- Six methods of prediction:
  - Compound covariate predictor (2 classes only)
  - Bayesian Compound covariate predictor (2 classes only)
  - K-nearest neighbor (2 or more classes)
  - Nearest centroid (2 or more classes)
  - Support vector machines (2 classes only)
  - Diagonal linear discriminant analysis (2 or more classes)


•

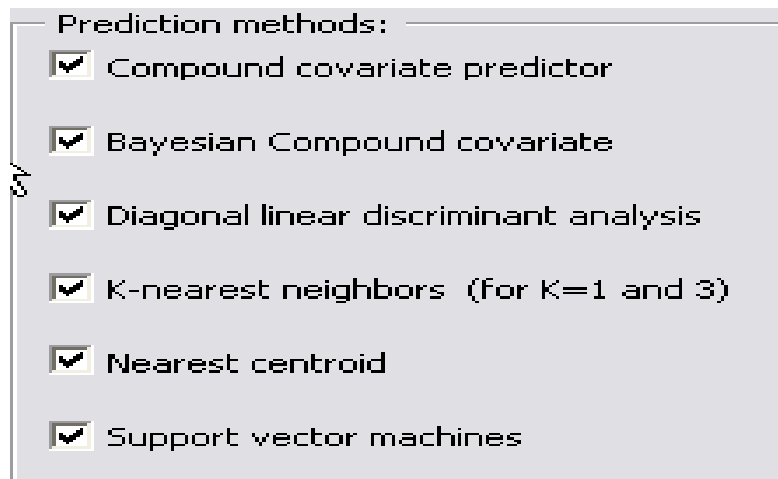



# C3. Cross-validating the classifier

- Leave-One-Out cross validation.
- K-Fold cross validation.
- $+0.632$  bootstrap cross-validation.
- Use leave-one-out cross-validation to compute a misclassification rate
- Re-compute the classifier, based on all but one sample
- Use the classifier to classify the sample which has been left out

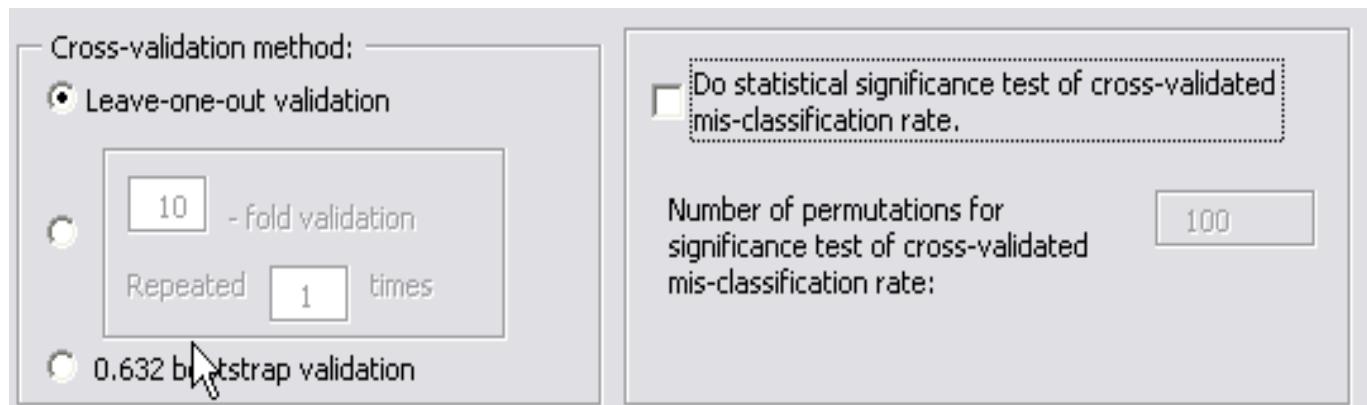

Cross-validation method:

- ☒ Leave-one-out validation
- ☐ 10 - fold validation  
Repeated 1 times
- ☐ 0.632 bootstrap validation

☐ Do statistical significance test of cross-validated mis-classification rate.

Number of permutations for significance test of cross-validated mis-classification rate: 100

# Permutation test

- Use a permutation test to assess the significance of the misclassification rate and univariate significance of each gene
- For each permutation of the class labels, re-run the cross-validation and obtain a new cross-validated misclassification rate
- The permutation p-value is based upon the rank of the misclassification rate using the original data, compared to all permutations

# Compound covariate predictor

- May only be used for classifying among two class labels
- Select genes which univariately classify the samples
- Form a compound covariate predictor as:

$$\sum_i t_i x_i \quad \left\{ \begin{array}{l} \text{where } t_i = \text{t-statistic, } x_i = \text{log-ratio,} \\ \text{and sum is taken over all significant genes} \end{array} \right.$$

- Determine the cutpoint of the predictor as the midpoint between its mean in one class and its mean in the other class

# Linear classifiers for two classes

$$l(\underline{x}) = \sum_{x \in F} w_i x_i$$

$\underline{x}$  = vector of log ratios or log signals

$F$  = features (genes) included in model

$w_i$  = weight for  $i$  - th feature

decision boundary  $l(\underline{x}) >$  or  $<$  cutoff

# Linear classifiers for two classes

- Diagonal linear discriminant analysis (DLDA)
- Compound covariate predictor
  - Bayesian compound covariate
- Support vector machine

# Diagonal linear discriminant analysis

- May be used for classifying among two or more class labels
- Use F-test to screen for genes which are univariately significant in classifying the samples
- Seeks a linear combination of the variables which has a maximal ratio of the separation of the class means to the within-class variance, where genes are assumed to be uncorrelated

# Bayesian Compound Covariate

- Compound Covariate score is computed for all the samples in the cross-validated training set.
- The CCP-scores of samples in each class of the training set are assumed to be from a Gaussian distribution.
- If prior probabilities are  $\frac{1}{2}$  - the BCCP is similar to the CCP.

# K-nearest neighbor

- May be used for classifying among two or more class labels
- Use F-test to screen for genes which are univariately significant in classifying the samples
- For  $k=1$  and  $k=3$ , finds the  $k$ -nearest neighbors in terms of Euclidean distance over only those genes which were univariately significant
- Classify based on the majority vote of the class labels of the  $k$ -nearest neighbors



# Nearest centroid

- May be used for classifying among two or more class labels
- Use F-test to screen for genes which are univariately significant in classifying the samples
- Compute the centroid of each class as a mean over all the training samples with that class label
- Classify test sample to be same class label as the nearest centroid, using Euclidean distance over only those genes which were univariately significant

# Support vector machines

(V. Vapnik)

- Implemented only for classifying among two class labels
- Select genes which univariately classify the samples
- The SVM predictor is implemented as a linear function of the log-ratios or the log-intensities over the significant genes, that best separates the data subject to penalty costs on the number of specimens misclassified.

# Class prediction tool

## Class prediction vs. binary tree prediction

- The class prediction tool has more options: may select all prediction methods simultaneously, may use paired samples, may use randomized variance option.
- The binary tree prediction tool splits the classes into groups of subclasses. At each node in the tree, the binary tree prediction tool decides how to split the classes into two groups based on either a leave-one-out or a K-fold cross-validation. The binary tree prediction tool may be useful if there is a hierarchical structure to the classes.
- However, the binary tree prediction may be very slow for a large number of samples. Therefore, a K-fold cross-validation should be used if the number of samples is large.
- Currently the tool is limited to five classes, and requires at least four samples per class for good prediction.

# Prediction Analysis Microarray

## PAM

- Uses Shrunk Centroid algorithm developed by Tibshirani's group (Stanford).
- Similar to Nearest Centroid but the centroids are shrunk towards each other based on shrinking the class means for each gene towards an overall mean.
- Amount of shrinking is determined by a tuning parameter  $\delta$  and the number of genes included in the classifier is determined by the value of  $\delta$ .

# Important notes

- Cross validation is only valid if the test set is not used in any way in the development of the model.
- With proper CV, the model must be developed from scratch for each leave-one-out training set. This means that feature selection must be repeated for each leave-one-out.

# [Hands-on instructions ]


[Class prediction –Pomeroy data]

1. Run the Class Prediction tool by clicking on **ArrayTools** → **Class prediction** → **Class prediction**.
2. Select the **Medulo vs Glio** variable as the column defining the classes. Check the box for using the Random Variance Model.
3. Choose the univariate significance  $\alpha=0.001$ .
4. Leave all other options at default levels, and click **OK**.
5. Note the Array Ids which have been misclassified by all methods.

# Part IV:

## RNA-Seq Count Data Importer

# RNA-Seq Count Data Importer

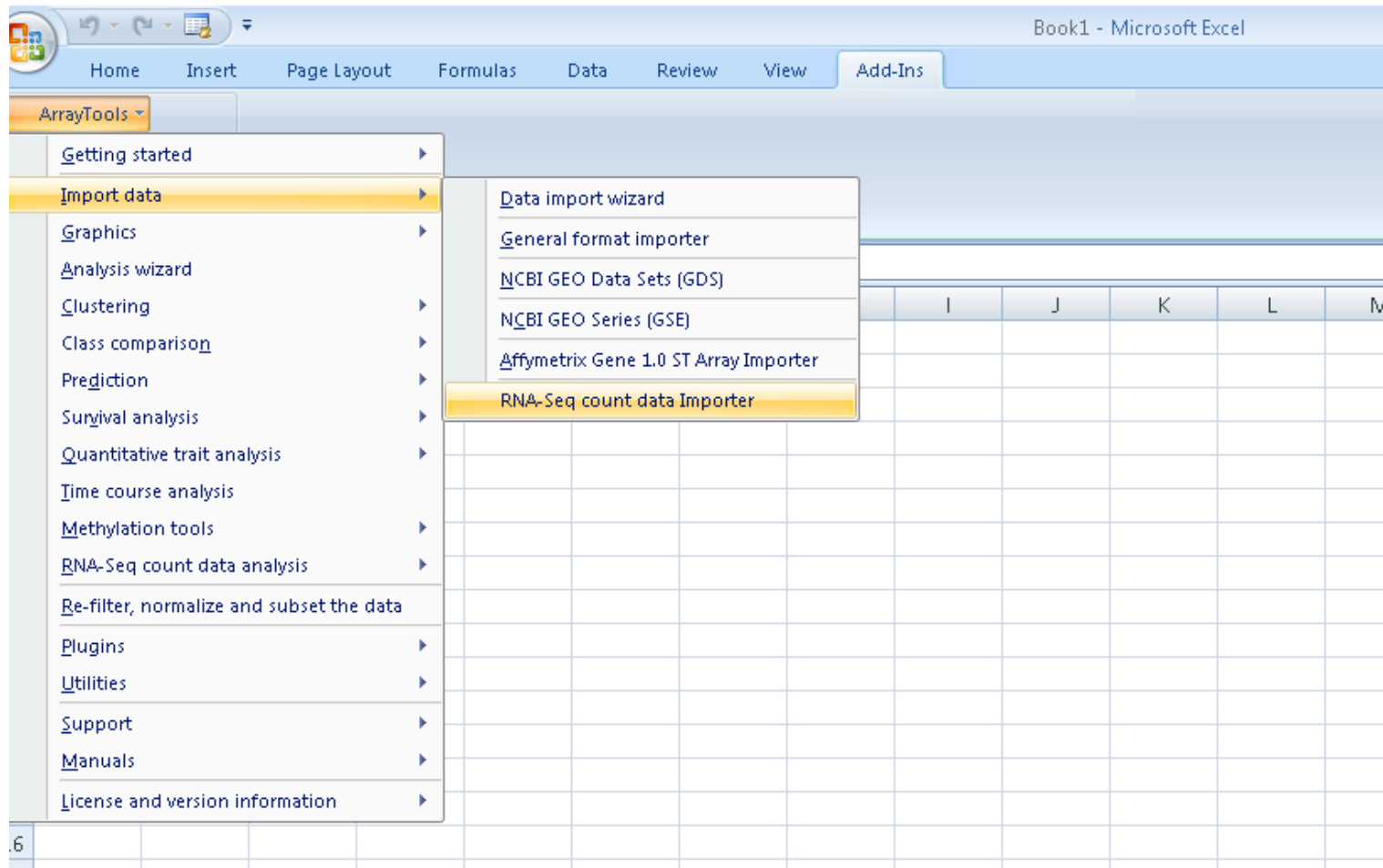



work - 10.0.6.181 - Remote Desktop Connection

CGHTools ▾  
VBAObjectManager ▾  
Menu Commands

Slides Outline

RNA-Seq Count Data Impoter

Select data file type

- ☒ I have count data files from HT-Seq
- ☐ Count data in one tab-delimited .txt file

Browse for data folder:

E:\data\GSE51403countable

Browse

Next Cancel

Sheet1 Sheet2 Sheet3

Ready

EN 100% 4:05 PM 5/3/2015

	A	B	C	D	E	F	G	H	I	J
1	Experiment	class								
2	Ctrl_Rep1	Ctrl								
3	Ctrl_Rep2	Ctrl								
4	Ctrl_Rep3	Ctrl								
5	Ctrl_Rep4	Ctrl								
6	Ctrl_Rep5	Ctrl								
7	Ctrl_Rep6	Ctrl								
8	Ctrl_Rep7	Ctrl								
9	E2_Rep1	E2								
10	E2_Rep2	E2								
11	E2_Rep3	E2								
12	E2_Rep4	E2								
13	E2_Rep5	E2								
14	E2_Rep6	E2								
15	E2_Rep7	E2								
16										
17										
18										
19										
20										
21										
22										
23										
24										
25										

Display the data		Ctrl_Rep1	Ctrl_Rep2	Ctrl_Rep3	Ctrl_Rep4	Ctrl_Rep5	Ctrl_Rep6	Ctrl_Rep7	E2_Rep1	E2_Rep2	E2_Rep3	E2_Rep4	E2_Rep5	E2_Rep6	E2_Rep7	Missi
1	UniqueID	1	2	3	4	5	6	7								
7	ENSG00000000938	1.191361	0.799476	0.982337	0.950698	0.959305	1.716867	1.78404	0	0	1.70454	1.233994	1.08408	0	1.616268	
8	ENSG00000000971	1.191361	0.799476	1.561364	1.518845	0.959305	0.817546	0.574129	0	0	0	0	1.08408	0	0	
19	ENSG000000001630	4.245764	6.312696	5.287283	5.072461	5.726299	5.37182	5.271872	5.726097	4.046669	6.005826	4.411602	5.288387	4.607452	5.403138	
23	ENSG000000002330	5.561097	7.004059	6.407813	6.34379	6.632482	6.432784	6.271453	6.712404	5.267583	6.665856	6.012393	6.390284	6.288069	6.33152	
26	ENSG000000002587	3.121361	3.597773	3.289978	2.912584	2.928011	3.722301	3.493692	2.931636	1.197762	2.464276	2.956218	2.722491	1.738857	2.623989	
28	ENSG000000002745	0	1.68769	0	0.950698	0	0.817546	0.983719	1.073792	2.287327	1.70454	2.68003	2.4542	1.116932	1.616268	
29	ENSG000000002746	2.278295	1.310885	0.982337	0	0.959305	0.817546	0	1.073792	1.197762	1.090599	0	1.696014	1.738857	1.023543	
36	ENSG000000003137	3.494344	3.409539	3.138374	3.493249	2.516484	3.230582	3.058607	4.135402	3.330107	4.253535	2.68003	3.608836	2.172026	4.043832	
42	ENSG000000003436	11.29676	10.62339	10.86383	10.87841	10.97749	10.93355	10.96868	9.713668	9.993291	9.806811	10.29906	9.752042	10.19643	10.1328	
52	ENSG000000004468	0	0	0.982337	0	0	0	0.983719	0	0	0	1.233994	0	0	0	
53	ENSG000000004478	12.37423	12.40696	12.56205	12.40151	12.50287	12.49133	12.55398	13.71324	13.45044	13.42375	13.67693	13.76099	13.57117	13.65068	
56	ENSG000000004660	6.237515	6.257912	6.052448	6.129856	5.896056	6.354294	6.097009	7.658943	7.275813	7.988889	7.25213	7.883894	7.223464	7.472705	
59	ENSG000000004776	2.278295	2.23338	2.555403	1.518845	1.938518	0.817546	0.983719	1.073792	0	1.70454	1.233994	2.722491	1.116932	1.023543	
60	ENSG000000004777	6.21276	7.73819	6.504002	7.157719	6.488177	7.037508	6.196513	7.864442	6.983018	8.432691	6.538348	7.73257	6.41206	7.000533	
62	ENSG000000004799	11.02276	10.76456	10.85474	10.90322	10.97207	10.93242	10.94097	9.803819	10.11779	9.877817	9.916553	9.651153	10.53781	10.00956	
64	ENSG000000004838	4.245764	5.433184	4.823677	4.075339	5.205007	4.622613	4.523251	4.99724	4.602234	4.924586	5.176183	5.586193	5.07587	4.745378	
66	ENSG000000004848	5.103317	4.582226	4.425522	4.29691	4.313786	4.622613	4.801304	4.893795	3.504708	4.071655	4.500469	4.403296	3.792583	4.292214	
70	ENSG000000004939	2.891109	3.305332	1.973424	2.721777	1.938518	2.266645	0.983719	2.109319	0	2.133851	2.956218	1.08408	2.172026	2.035172	
71	ENSG000000004948	5.350288	5.618546	5.13538	5.187717	4.920081	4.530431	5.196954	7.794362	8.183979	7.658734	7.892575	7.950731	7.998929	7.994888	
74	ENSG000000005001	8.120367	9.293325	8.647538	9.018781	8.491298	9.039735	8.4632	8.308421	7.770595	8.53289	7.26487	8.099535	7.683341	7.668526	
83	ENSG000000005108	1.834857	1.310885	0.982337	1.518845	0.959305	0	0	1.073792	0	1.70454	0	0	1.116932	1.023543	
89	ENSG000000005206	7.866328	9.37443	8.369703	8.674677	8.348126	8.75917	8.24813	8.992318	8.042474	9.508742	7.772809	8.849653	7.972317	8.149828	
95	ENSG000000005379	1.834857	1.310885	0.982337	1.518845	0.959305	0	0	1.073792	0	1.70454	0	0	1.116932	1.023543	

Experiment descriptors

Gene annotations

Filtered log intensity

Gene identifiers

	A	B	C	D	E	F	G	H
	UniqueID (Double-click)	Name	Accession	UGCluster	Symbol	EntrezID	Chromosome	Cytoband
7	<a href="#">ENSG00000100000</a>	FGR proto	NM_001044	Hs.1422	FGR	2268	1	1p36.2-p3
8	<a href="#">ENSG00000100000</a>	complement	NM_000168	Hs.363396	CFH	3075	1	1q32
19	<a href="#">ENSG00000100000</a>	cytochrome	NM_000788	Hs.417077	CYP51A1	1595	7	7q21.2
23	<a href="#">ENSG00000100000</a>	BCL2-asso	NM_004322	Hs.370254	BAD	572	11	11q13.1
26	<a href="#">ENSG00000100000</a>	heparan s	NM_005111	Hs.507348	HS3ST1	9957	4	4p16
28	<a href="#">ENSG00000100000</a>	wingless-1	NM_016068	Hs.272375	WNT16	51384	7	7q31
29	<a href="#">ENSG00000100000</a>	HECT, C2 s	NM_001288	Hs.164453	HECW1	23072	7	7p13
36	<a href="#">ENSG00000100000</a>	cytochrome	NM_001275	Hs.91546	CYP26B1	56603	2	2p13.2
42	<a href="#">ENSG00000100000</a>	tissue fact	NM_001035	Hs.516578	TFPI	7035	2	2q32
52	<a href="#">ENSG00000100000</a>	CD38 mole	NM_001775	Hs.479214	CD38	952	4	4p15
53	<a href="#">ENSG00000100000</a>	FK506 bin	NM_002001	Hs.524183	FKBP4	2288	12	12p13.33
56	<a href="#">ENSG00000100000</a>	calcium/c	NM_032295	Hs.8417	CAMKK1	84254	17	17p13.2
59	<a href="#">ENSG00000100000</a>	heat shock	NM_144611	Hs.534538	HSPB6	126393	19	19q13.12
60	<a href="#">ENSG00000100000</a>	Rho GTPas	NM_001117	Hs.515364	ARHGAP3	115703	19	19q13.12
62	<a href="#">ENSG00000100000</a>	pyruvate c	NM_002611	Hs.8364	PK4	5166	7	7q21.3
64	<a href="#">ENSG00000100000</a>	zinc finger	NM_015895	Hs.526735	ZMYND10	51364	3	3p21.3
66	<a href="#">ENSG00000100000</a>	aristaless	NM_139055	Hs.300304	ARX	170302	X	Xp21.3
70	<a href="#">ENSG00000100000</a>	solute car	NM_000344	Hs.443948	SLC4A1	6521	17	17q21.31
71	<a href="#">ENSG00000100000</a>	calcitonin	NM_001116	Hs.489127	CALCR	799	7	7q21.3
74	<a href="#">ENSG00000100000</a>	protease,	NM_022111	Hs.459709	PRSS22	64063	16	16p13.3
83	<a href="#">ENSG00000100000</a>	thrombos	NM_015201	Hs.120855	THSD7A	221981	7	7p21.3
89	<a href="#">ENSG00000100000</a>	signal per	NM_001015	Hs.744026	SPPI2R	56928	19	19p13.3
<div> <div>Experiment descriptors</div> <div>Gene annotations</div> <div>Filtered log intensity</div> <div>Gene id</div> </div>								
Ready Filter Mode								

# Differential expression analysis of RNA-Seq Count Data

The screenshot shows the ArrayTools software interface. The 'ArrayTools' menu is open, highlighting 'RNA-Seq count data analysis'. Below this menu, a data table is displayed with columns E through N. The table contains gene expression data, including gene symbols, Entrez IDs, chromosome numbers, cytobands, GO terms, and Ensembl IDs. A 'Filter' column is highlighted in yellow, and a 'TRUE' value is visible in the first row of the data table.

	E	F	G	H	I	J	K	L	M	N
							Filter			
	Symbol	Entrez	Chrom	Cytoba	GO	Ensembl				
	FGR	2268	1	1p36.2-p3	#####	ENSG0000	TRUE			
	96 CFH	3075	1	1q32	#####	ENSG0000	TRUE			
	77 CYP51A1	1595	7	7q21.2	#####	ENSG0000	TRUE			
					#####	ENSG0000	TRUE			
					#####	ENSG0000	TRUE			
					#####	ENSG0000	TRUE			
	53 HECW1	23072	7	7p13	Biological	ENSG0000	TRUE			
	5 CYP26B1	56603	2	2p13.2	#####	ENSG0000	TRUE			
	78 TFPI	7035	2	2q32	#####	ENSG0000	TRUE			
	14 CD38	952	4	4p15	#####	ENSG0000	TRUE			
	33 FKBP4	2288	12	12p13.33	#####	ENSG0000	TRUE			
	59 ENSG0000	84254	17	17p13.2	#####	ENSG0000	TRUE			
	59 ENSG0000	126393	19	19q13.12	Biological	ENSG0000	TRUE			
	50 ENSG0000	115703	19	19q13.12	#####	ENSG0000	TRUE			
	50 ENSG0000				#####	ENSG0000	TRUE			

# Find DE genes with edgeR

The screenshot shows a web-based interface titled "Differential gene expression with edgeR". It contains several input fields and dropdown menus for configuring the analysis. The "Class variable column" is set to "class". The "Find differentially expressed genes using" method is set to "Classic edgeR approach". The "Find differentially expressed genes by" criterion is set to "threshold p-values". The "Threshold p-value for significance" is set to "0.001". The "Threshold fdr for significance" is set to "0.05". The "Output name" is set to "EDGER". There are two buttons at the bottom: "Submit" and "View README".

Class variable column	class
Find differentially expressed genes using	Classic edgeR approach
Find differentially expressed genes by	threshold p-values
Threshold p-value for significance	0.001
Threshold fdr for significance	0.05
Output name	EDGER

**Submit**

[View README](#)

# Output Results File

variable: class  
Ctrl  
method: Classic edgeR  
termination of significance by: p-values  
number of genes: 62069

ult:

number of identified genes: 5184

UniqueID	Anno- tations	Symbol	Name	EntrezID	Accession	UGCluster	Ensembl	DefinedGeneList	P_Value	FDR	log2FC
ENSG00000196208	<a href="#">Info</a>	<a href="#">GREB1</a>	growth regulation by estrogen in breast cancer 1	<a href="#">9687</a>	<a href="#">NM_014668</a>	<a href="#">Hs.467733</a>	<a href="#">ENSG00000196208</a>	Downregulated of MTA-3 in ER-negative Breast Tumors	< 1e-4	< 1e-4	4.2848
ENSG00000147256	<a href="#">Info</a>	<a href="#">ARHGAP36</a>	Rho GTPase activating protein 36	<a href="#">158763</a>	<a href="#">NM_001282607</a>	<a href="#">Hs.22905</a>	<a href="#">ENSG00000147256</a>		< 1e-4	< 1e-4	6.041
ENSG00000107562	<a href="#">Info</a>	<a href="#">CXCL12</a>	chemokine (C-X-C motif) ligand 12	<a href="#">6387</a>	<a href="#">NM_000609</a>	<a href="#">Hs.522891</a>	<a href="#">ENSG00000107562</a>	CXCR4 Signaling Pathway, Pertussis toxin-insensitive CCR5 Signaling in Macrophage, Axon guidance, Chemokine signaling pathway, Cytokine-cytokine receptor interaction, Intestinal immune network for IgA production, Leukocyte transendothelial migration, Rheumatoid arthritis	< 1e-4	< 1e-4	3.5643

www.ncbi.nlm.nih.gov/sites/entrez?db=gene&cmd=search&term=9687

# Part V: Plug-ins



# Plug-in utility

- A plug-in utility now allows users to create their own tools by writing their own scripts written in the R language
- Tools created using the plug-in utility can be distributed to other users, and added to the Plugin menu
- The user-created plug-ins are stored in the Plugins folder of the ArrayTools installation folder

# Included plugins

- Analysis of Variance – Up to four-way ANOVA. Options to include blocking factors or use random variance model.
- ANOVA of log intensities – For dual-channel non-reference designs, model includes gene-specific array effect, dye effect, and class effect. Option to use random variance model.
- ANOVA for Mixed Effects Model – Allows up to three fixed effects and one random effect.
- M vs A plot – For dual-channel data, plots log-ratio vs average log-intensity for all arrays.
- Pairwise correlation – Plots heat map showing the matrix of pairwise correlations among all arrays.
- Smoothed CDF – Plots smoothed cumulative distribution function of log-red and log-green, or log-ratio for all arrays.
- Export 1- and 2-color data to R – Exports data from Project Workbook to files which can be imported into R.

# [Additional Plugins]

- Class Prediction using TopScoring Pairs: This plugin is a different tool for class prediction by using the top-scoring pairs (TSP) classifier developed by Geman et al.
- Random Forest: This tool is another alternative to class prediction and the random forest is built from the ensemble learning method - methods that generate many classifiers and aggregate their results. The random forest is robust against overfitting and has been demonstrated to have performance competitive with the other classifiers.
- TimeSeries: This plug-in can be used for regression analysis of time series expression data.

**Create Plug In**

FileNames

R-Script Full Path:  Browse...

Plug In Filename:

Plug In Title:

Plug In Description:

Data to Send to R-Script

☐ Either Filtered Normalized Log Intensity or Filtered Normalized Log Ratio

☐ Experiment Design Worksheet

☐ Gene Identifiers Worksheet

Variable Names

☐ Two Color Unnormalized Intensities

☐ Filtered Normalized Log Ratio

☐ One Color Unnormalized Log Intensity

☐ Filtered Normalized Log Intensity

Cancel Without Saving

Create Plug In

## Part VI:

Independent practice  
(if time permits)

# Further help

- We hope this class has been helpful to you. This class was not designed to be comprehensive, but only an introductory overview of the features in BRB-ArrayTools. More information about the software may be obtained from the User's Manual (may be viewed by clicking on **ArrayTools -> Support -> Manuals -> User's Manual**).
- Supplementary material on analysis algorithms may be found in the BRB technical reports:  
<http://linus.nci.nih.gov/~brb/TechReport.htm>

# Acknowledgements

- Dr. Richard Simon and Biometrics Research Branch staff members.
- BRB-ArrayTools development team (past and present).
- User community!!

# Technical support

- For questions of a general nature, post a message to the BRB-ArrayTools Message Board:  
<http://board.emmes.com/phpBB3/index.php>
- To report bugs, send email to [arraytools@emmes.com](mailto:arraytools@emmes.com)



# Feedback on this class

- Please fill out a feedback form before you leave the class.
- Please make your comments specific enough to enable us to adjust this presentation for future classes.
- Thank you for participating in this class!!

# Exercise Section

**U**sing the breast tumors sample data set, find genes that are differentially expressed for patients before and after treatment:

- Obtain a gene list that contain no more than 40% of False discoveries with 95% confidence.
- Choosing an alternative method to the Multivariate Permutation test to control for false discoveries obtain another gene list with a 95% confidence level and controlling for 40% False discoveries.
- Using all genes in this sample dataset, run a scatter plot of phenotype averages with 1.5 fold difference and comment on the up/downward regulated genes.

**BRB-Digital Gene Expression (BRB-DGE)**


***<http://linus.nci.nih.gov/bdge/>***



# BRB-Digital Gene Expression (BRB-DGE)

- An RNA-Seq pre-processing tool to produce count data file
- Easy accessible and user-friendly for general scientists
- Includes well-known software tools for RNA-Seq data pre-processing

# BRB-DGE Workflow

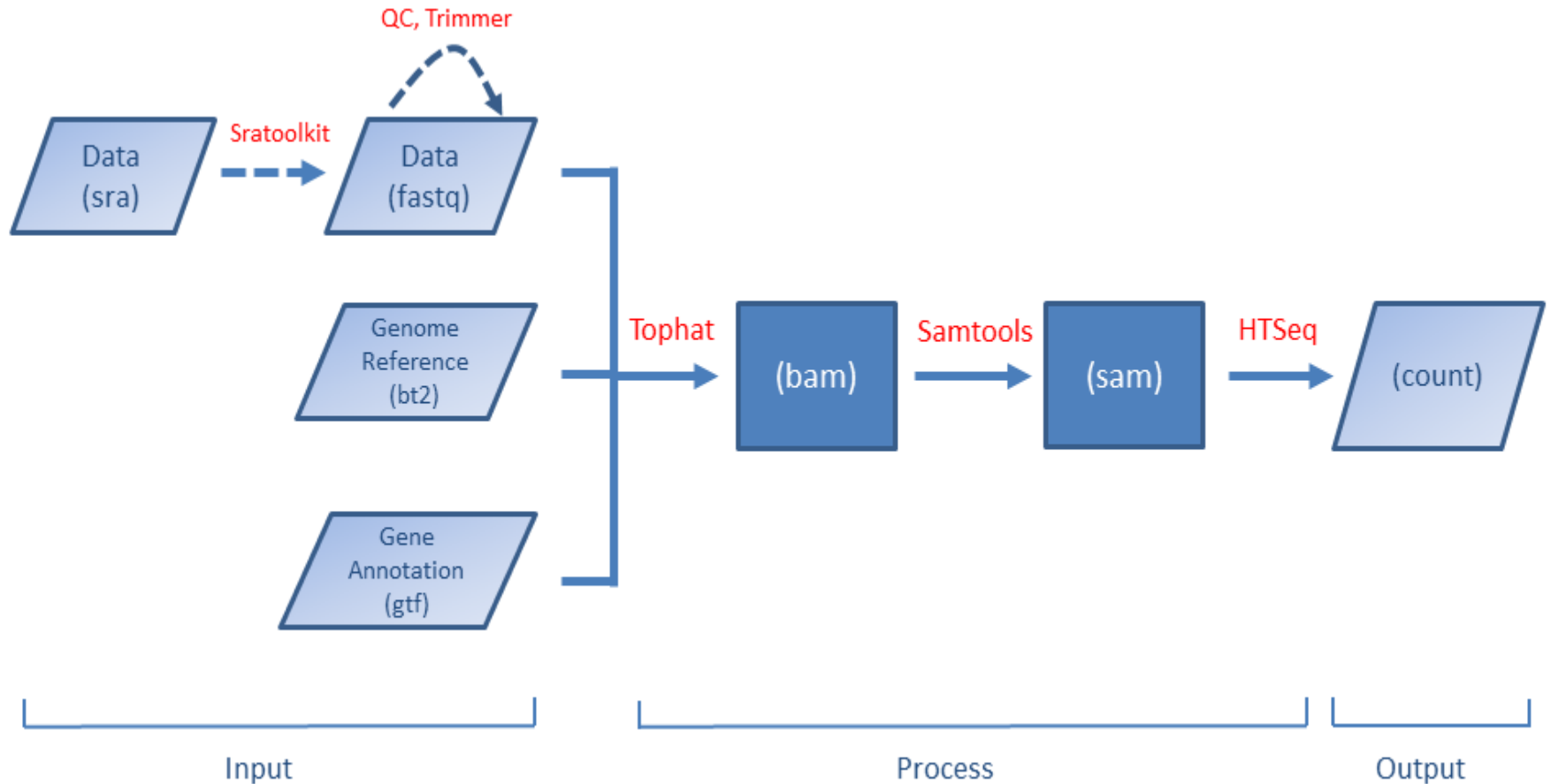


# OS and Hardware Requirements

- Ubuntu 64-bit Desktop OS (version 12.04 or 14.04). Can be downloaded at <http://www.ubuntu.com/download/desktop>
- For non-ubuntu OS computers, virtual machine or dual boot can be used in order to install Ubuntu OS

# Required Input Files

- fastq - sequence data file
- Reference genome files in bt2 format
- Reference genome file in fa format (optional but recommended)
- Gene annotation file in gtf format
- samples.txt

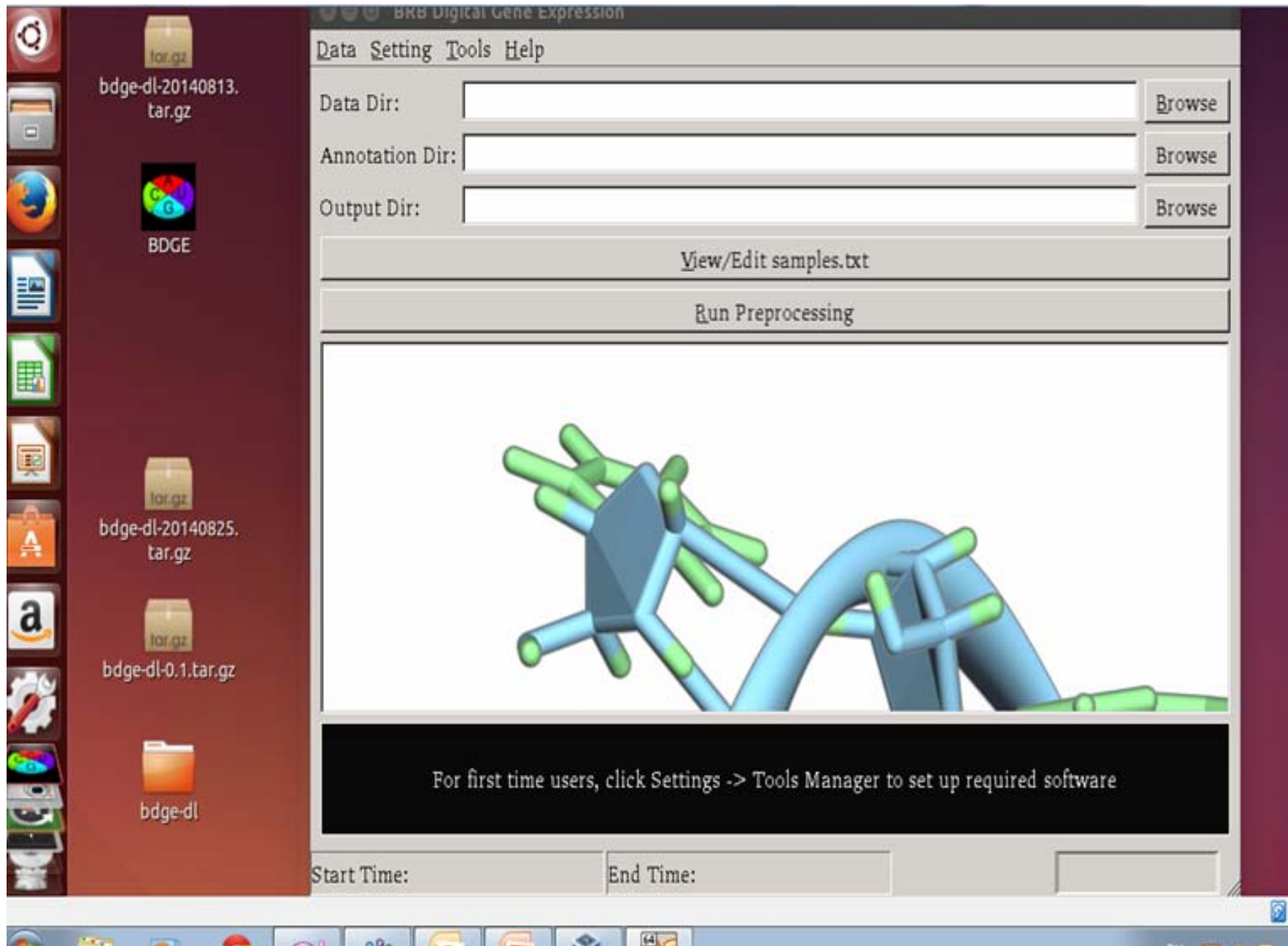



# Output Count File

A1										
	A	B	C	D	E	F	G	H	I	J
1	A1BG	16								
2	A1BG-AS1	12								
3	A1CF	1								
4	A2M	1213								
5	A2M-AS1	2								
6	A2ML1	6								
7	A2MP1	1								
8	A3GALT2	0								
9	A4GALT	129								
10	A4GNT	0								
11	AA06	0								
12	AAAS	965								
13	AACS	1017								
14	AACSP1	0								
15	AADAC	0								
16	AADACL2	0								
17	AADACL2-	0								
18	AADACL3	0								
19	AADACL4	0								